MULTITHERMO

Mode d'emploi - Instruction Manual
Gebrauchsanweisung - Instrucciones de uso
Manual de instruções

Terraillon®

Français

THERMOMETRE INFRAROUGE 3 EN 1

Cher(e) client(e),

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3-IN-1 INFRARED THERMOMETER

Dear Customer,

Thank you for buying this Terrailon infrared electronic thermometer. We hope it gives you excellent service and to ensure complete satisfaction, we recommend you read these instructions carefully.

The infrared sensor of the Multithermo thermometer records the heat generated by the eardrum (this reflects the temperature of the nervous system in the brain which regulates the temperature of the body) or on the surface of the forehead. It also measures the ambient temperature of your room.

Quick and accurate, it takes 512 measurements per second and displays the highest temperature in 1 second. Practical and hygienic, the sensor is washable so a protective tip does not need to be used.


It's easy to use, switching from forehead measurement mode to ear measurement mode when you remove the protective cap and its liquid crystal display with big figures is easy to read. It's discreet but always there within reach on its stand, displaying the time and ambient temperature on standby.


IMPORTANT

1. Temperature varies according to the part of the body where it is taken. That's why temperature readings taken orally, under the arm, rectally, in the ear or on the forehead should not be compared.
2. Temperature also varies from one person to another, according to age and the time of day.

TIP

Take your temperatures when you are completely well. You can use these results as reference temperatures and compare them with your temperature readings when you are ill. Is your temperature abnormally high or abnormally low compared to normal? Do not hesitate to consult a doctor if in doubt.

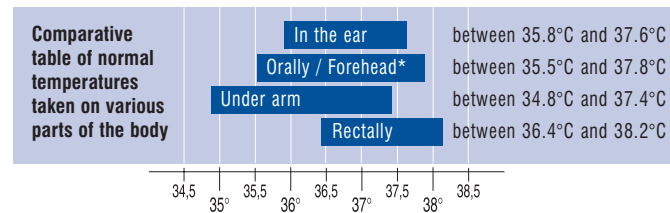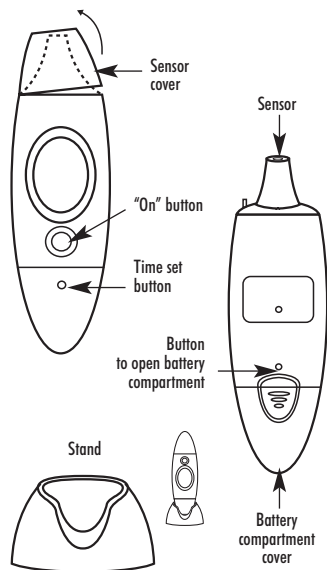USING YOUR THERMOMETER FOR THE FIRST TIME

The device is supplied with a 3V battery (CR2032) already in place, so it is ready to use.

SETTING THE TIME

The first time you use this thermometer or after changing the battery, you have to set the time.

Open the battery compartment using the tip of a biro. Press the time set button. The hour numbers flash. Press the 'on' button of the thermometer to set the

hour. Confirm by pressing the time set button. The minute numbers flash. Press the 'on' button of the thermometer to set the minutes. Confirm by pressing the time set button.

On standby, the screen will alternately display the ambient temperature and the time.

TAKING TEMPERATURES

With its protective cap on, the device is in forehead mode by default. When the protective cap is removed, the thermometer automatically switches to ear mode.

1. Forehead

Press the temperature reading button. There is a beep, then the latest temperature in the memory is displayed. Press the temperature reading button and keep it depressed. Place the sensor on one of your temples. Slide the sensor across your forehead to the other temple. Release the button to stop the measurement. There is a beep and your temperature is displayed.

**The forehead temperature shown corresponds exactly to the oral temperature.*



The forehead temperature can be subject to the influence of outside climatic conditions. In case of doubt, do not hesitate to take the temperature in the ear.

2. Ear

Remove the protective cap. The device automatically switches to ear mode. Press the temperature reading button.


There is a beep indicating that the device is ready, then the latest temperature in the memory appears. Press the temperature reading button and keep it depressed. Insert the sensor in the auditory canal until a beep informs you that the temperature reading is finished. Release the button. Remove the sensor from the ear. Your temperature is displayed.

Important: In order to get an accurate measurement, make sure that you position the sensor in the ear correctly so that the tip of the sensor picks up the heat generated by the eardrum. To do this, you can pull the ear to straighten the canal and turn the sensor slightly inside (For users under the age of 1, pull the ear backwards gently; for children over 1, pull the ear up and back).

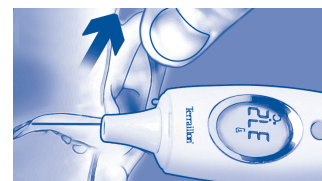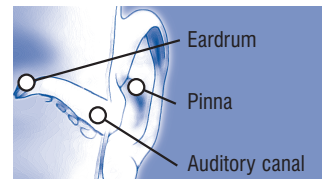CHANGING THE BATTERY

When the "  " symbol appears on the screen, this means that the battery is running low. When the screen displays "Lo", this means that the battery has run out and has to be changed.

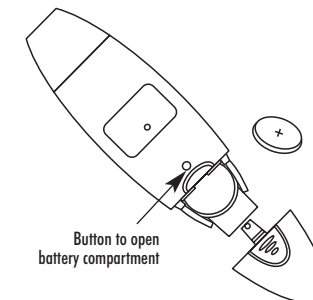Open the battery compartment by pressing the button on the back of the device with the tip of a pen. Remove the used battery.

Replace it with a 3V lithium CR2032 battery, inserting it into the designated space.

The widest part of the battery must be on top. The "+" must remain visible once the battery is in place.

Replace the cover. The device is ready to use.

(Do not throw away used batteries. Put them in a special container for collection and recycling). 

SPECIAL MESSAGES

	Temperature reading in progress
	Ear temperature reading mode
	Forehead temperature reading mode
	Clock
	Temperature is over 50°C
	Temperature is over 10°C
	Memory. Result of latest temperature reading.

 This symbol appears when the battery is starting to run low. You can still use the appliance, but when the "Lo" sign appears, the battery must be changed because the device can no longer function normally.

 Degrees Celsius

 Temperature exceeds operating conditions (16-40°C)

SWITCHING OFF YOUR THERMOMETER

The thermometer goes on standby automatically after 30 seconds if no buttons are activated.

LOOKING AFTER YOUR THERMOMETER

Clean the thermometer sensor with cotton wool soaked in surgical spirit before and after each use. Do not use abrasive substances. Do not immerse the thermometer in water or any liquid whatsoever because only the sensor is impermeable.

WARNINGS

1. Taking your temperature is a way of monitoring it. It does not mean a diagnosis can be made or a treatment prescribed. A doctor should always be informed of abnormal readings. Do not change the doses of drugs prescribed by the doctor under any circumstances.
2. Do not squeeze, drop or break the thermometer.
3. Keep it out of reach of children.
4. Do not expose the thermometer to the sun and do not store it close to a source of heat but keep it somewhere clean and dry.
5. Do not immerse the thermometer in water or any liquid whatsoever because only the sensor is impermeable.
6. The Multithermo thermometer can only be used to take the temperature in the ear and on the forehead and to measure ambient temperature. Do not use it to take the temperature rectally, orally or under the arm.
7. Do not touch the end of the sensor with the fingers.
8. Do not take a temperature reading less than 30 minutes after eating, exercising or taking a bath.
9. Take temperature readings at one-minute intervals. Two successive temperature readings may give slightly different results. In this case, use the average of the two.

MANUFACTURER

AVITA CORPORATION
9F, No.78, Sec.1,
Kwang-Fu Road,
San-Chung, 241
Taipei County
Taiwan, R.O.C
CE0197

TECHNICAL DATA

Range of measurement :	10°C - 50°C
Tolerance :	Between 36 and 39°C : +/-0.2°C Under 36°C : +/- 0.3°C Over 39°C : +/- 0.3°C
Graduation :	0.1°C
Operating conditions :	16 to 40°C 15-95% RH
Storage conditions :	-20 to 50°C 15-95% RH
Battery :	3V CR2032

GUARANTEE

This product is guaranteed for 2 years against defects in materials and manufacturing defects. During this period, such defects will be repaired free of charge (proof of purchase required in the case of a claim under the guarantee). This guarantee does not cover damages resulting from accidents, improper use or negligence. In the case of a claim, your first point of contact should be the shop where you bought your thermometer.

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