MultiThermo 50

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FR	Thermomètre infrarouge 3 en 1	p 2
EN	3 in 1 infrared thermometer	p 7
DE	3-in-1-Infrarotthermometer	p 12
IT	Termometro a raggiinfrarossi 3 in 1.....	p 17
ES	Termómetro infrarrojo 3 en 1	p 22
NL	3-in-1-infraroodthermometer.....	p 27

THERMOMETRE INFRAROUGE 3 EN 1

Cher(e) client(e),

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3 IN 1 INFRARED THERMOMETER

Dear Customer,

Thank you for purchasing the Terraillon infrared electronic thermometer. We hope you will be entirely satisfied with your purchase. To ensure your satisfaction, please read the information leaflet thoroughly before using this product.

Multithermo 50's infrared sensor detects and measures the heat on the forehead or generated by the eardrum (it shares the same blood supply with the temperature control centre in the brain and accurately reflects the body's core temperatures). It can also measure the room's ambient temperature. With 512 scans per second, it is fast and accurate. The highest temperature reading is displayed after one second and up to 10 measurements can be saved under its memory feature. Practical and hygienic, its reusable lens is perfectly washable, no more disposable lens filters. Easy to use, it can be automatically adapted for measuring forehead or ear temperatures just by removing the protective cap. Its liquid crystal display with large numbers is easy to read. When not in use, it can be discretely stored on its base for easy access while displaying the time and room temperature.


READ BEFORE USING

Important


1. Temperature will vary depending on the part of the body. This is why you should not compare axillary (armpit), rectal, auricular (ear) or forehead temperature readings.
2. Temperature will also vary depending on the person, age, and time of day.

Recommendation

Determine your normal temperature range by taking your temperatures when healthy. These readings can be used as a form of reference to compare with future readings when ill. Is your temperature abnormally high or low in comparison to your normal temperature? Contact your physician immediately to be sure.

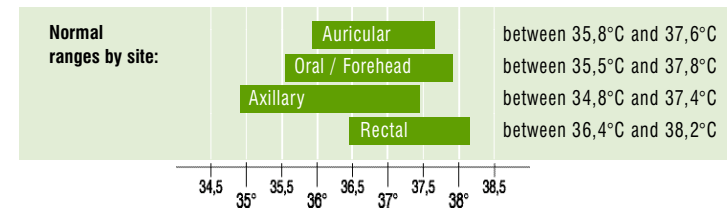PRODUCT DESCRIPTION

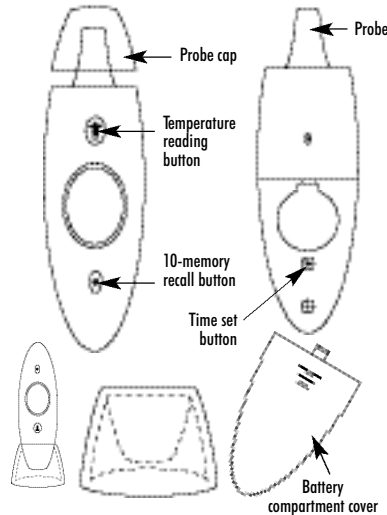FIRST USAGE

1. Power

A 3V (CR2032) battery is included inside unit. Unit is ready for immediate use. Please remove the battery tab before use.

2. Setting the Clock

Set the clock when using for the first time or after changing the battery.

- Open the battery case by sliding off the cover on the back of the unit.
- Push the "time setting" button. The hour numbers will start blinking.
- Push the "⊕" button to select.
- Push the "time setting" button again to set. The minute numbers will start blinking.

- Push the "⊕" button to select.
- Push the "time setting" button again to set. While not in use, the screen will display the time and room temperature.

3. How to choose the temperature measuring unit (F° or C°)? p. 11.

TAKING TEMPERATURE


When the protective cap is attached, the unit is automatically in forehead measure mode. It is easily adapted for the ears by removing the protective cap.

1. Forehead Mode

- The cap attaches by fitting over the probe.
- Press the "⊕" button. All screen symbols will appear (fig. 1), followed by the last temperature taken (fig. 2).
 - When you hear a beep, press and hold down the "⊕" button.
 - Hold the sensor against one of your temples.
 - Slide the sensor across your forehead to the other temple.
 - Release the button to stop measuring. You will hear a beep and the temperature will be displayed (fig. 3).

The displayed forehead temperature would be the same as the oral temperature.


Exterior conditions can influence forehead temperature. If you are unsure of the results, try taking your ear temperature.

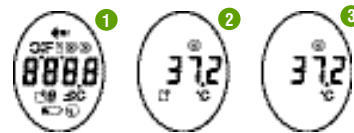
2. Ear Mode

Remove the protective cap. The unit is now in ear mode.

- Press the "⊕" button. All screen symbols will appear (fig. 1), followed by the last temperature taken (fig. 2).
- When you hear a beep the thermometer is ready (fig. 3), press and hold down the "⊕" button. (fig. 4)
- Fit the sensor probe into the ear canal until a beep signals the end of the measuring process.
- Release the button.
- Remove sensor from ear canal. Your temperature will be displayed.

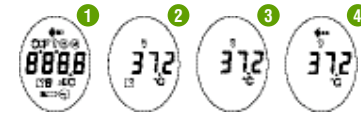
Important: To make an accurate reading, make sure that the sensor is securely positioned in the ear canal so that the sensor lens registers heat generated by the eardrum. To do this, gently pull your ear to straighten the ear canal and fit sensor snugly by turning lightly. (The thermometer should not be used as an ear thermometer on babies between the age of 0-12 months.)

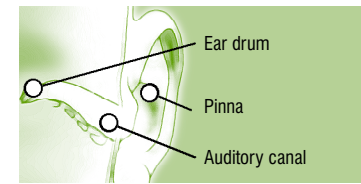MEMORY

To display the last temperature taken, press the "⊕" button, the temperature will appear along with its memory number.

To display the previous memories, press the "⊕" button and the memory numbers will appear.

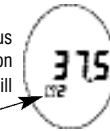CHANGING BATTERIES

When the symbol "⚡" appears, this means that the battery is low. When the display shows the symbol "Lo" as well and you hear a warning beep, this means that the battery no longer works and should be changed. Press any button to stop warning beep.

- Open the battery case by sliding off the cover on the back of the unit.
- Remove used battery.
- Replace with a lithium 3V CR2032 battery in designated area. The larger part of the battery should be facing up. You should still be able to see the "+" sign when battery is installed.
- Slide the cover back on. Unit is ready for immediate use.

Do not throw away used batteries. Take them to your local battery collection center to be recycled.

DISPLAY SYMBOLS

	Scanning
	Ear Mode
	Forehead Mode
H	Temperature higher than 50°C
Lo	Temperature lower than 10°C
	Memory. Last Temperature Reading
	Clock
C	Reading in Celsius

This symbol appears when the battery is low. The thermometer can still be used, but when you see the symbol "Lo", the battery should be changed because the thermometer can no longer function properly.

Err Ambient operating temperature (16-40°C) out of range

If the "👂" or "😊" is blinking and you hear a beeping signal, the thermometer is not ready for use. Push one of the buttons then wait for the symbol to stop blinking.

STAND BY MODE

The thermometer goes into standby automatically after 30 seconds of inactivity.

MAINTENANCE

Clean thermometer sensor with a cotton swab moistened with alcohol before and after each usage. Do not use abrasive cleaners. Do not submerge thermometer in water or any other liquid because only the sensor is waterproof.

PRECAUTIONS

1. Use of this thermometer is for monitoring purposes only and is not intended as a diagnostic tool. Make sure that any unusual readings are communicated to your physician. Under no circumstances should you change your prescribed medicine doses.
2. Do not grip strongly, drop, or crush thermometer.
3. Keep out of reach of children.
4. Keep thermometer away from direct sunlight and store well away from any source of heat in a clean, dry location.
5. Do not submerge thermometer in water or any other liquid because only the sensor is waterproof.
6. Multithermo 50 is only meant for taking temperatures in the ear, on the forehead and in the room. Do not use for taking rectal, oral, or axillary temperatures.
7. Do not touch the sensor lens with fingers.
8. Do not take temperature within 30 minutes of eating, exercising, or taking a bath.
9. Wait one minute between each temperature reading. Two readings right after each other can give slightly different measurements. Should this occur, use the average of the two.

How to choose the temperature measuring unit (Fahrenheit F° or Celsius C°)?

When the appliance is turned off: press button "⊕" to start it up. The F° and C° symbols will appear on the screen. Rapidly press button "⊕" (before the last recorded temperature comes onto the screen) to change from one measuring unit to the other. The thermometer is ready for use when the memory symbol "M" is displayed.

MANUFACTURER

AVITA CORPORATION
No.200.205 Provincial
Free way Song Ling Town,
Wu Jiang City,
Jiang Su, P.R.C. ZIP:215200
CE0197

TECHNICAL SPECIFICATIONS

Operating temperature range:	10°C - 50°C
Tolerance level:	between 36 and 39°C: +/- 0.2°C Below 36°C: +/- 0.3°C Above 39°C: +/- 0.3°C
Graduation:	0.1°C
Operating conditions:	16 to 40°C 15-95% RH
Storing conditions:	-20 to 50°C 15-95% RH
Battery:	1 x 3V CR2032

GUARANTEE

This unit is guaranteed for 2 years against any material or manufacturing defects. During this period, any such defects will be repaired free of charge (proof of purchase must be furnished for each claim). This guarantee does not cover damage due to accidents, improper use, or negligence. Should you like to make a claim, please first contact the store where you made your original purchase.

This Multithermo 50 probe is antibacterial and has been tested by an approved laboratory.

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HANSON (UK) LTD
TERRAILLON ASIA PACIFIC LTD