What is NAIS THERAPY?

NAIS THERAPY is a muscle and nerve stimulator which allows you to treat shoulder and neck pain, strains, bruises, sore muscles, etc. yourself. The effects on the muscles include thorough relaxation and gentle pain relief.

This device is subject to technical inspection by TÜV Rheinland Product Safety GmbH and bears the TÜV/GM mark.

NAIS THERAPY is based on gentle stimulation transferred to the hurting and tense body parts by way of "electronically controlled pads". This therapeutic stimulation penetrates deep into the tissue where it relaxes the muscles thoroughly and improves blood circulation, thereby usually also providing noticeable pain relief.

Simply place the NAIS pads on the parts of your body that you wish to treat. The pads are flat and convenient. A hand-held unit facilitates convenient control of the therapeutic stimulation. Power supply: 2 x 1.5 V alkaline batteries, type: LRO3/micro.

The NAIS pads

The NAIS THERAPY stimulator is supplied with one set of selfadhesive pads (2 in total). Note that the treatment always requires two pads and will not work with only one.

The NAIS pads will work more efficiently and last longer if you apply them correctly and make sure your skin is clean and free from grease. It is therefore in your own interest that you should clean the relevant parts of your skin and remove any body hair, if necessary, before using NAIS THERAPY.

1. Remove protective film

Before removing the protective film, connect the NAIS pad and the cable with the push contact, situated in the upper section of the pad. Then remove the protective film.

2. Connect NAIS pads to the hand-held control unit

When connecting the NAIS pads to the hand-held control unit, please note that the connector mechanism is designed in compliance with the 93/42/EEC directive on medical products, thus preventing unwanted contact with other sources of electrical power or electronic signals.

Never turn on the hand-held control unit before both pads are in place. Never stick the adhesive sides of two pads together!


3. Position the NAIS pads

See suggestions for treatment,

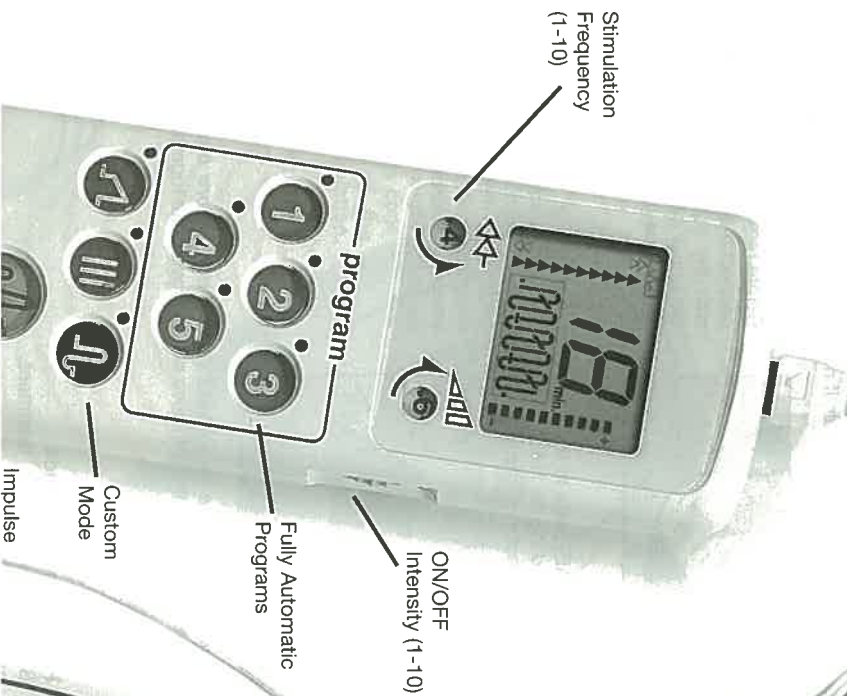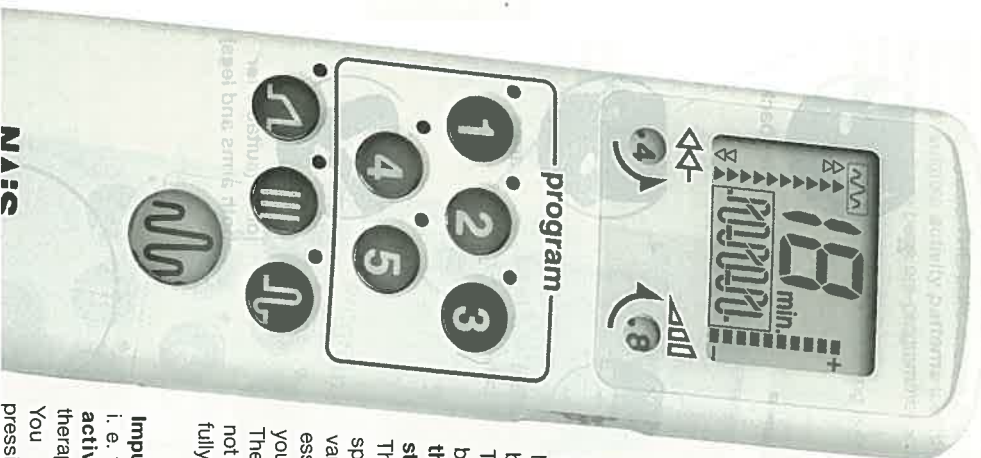
Simple electronic control of NAIS stimulation therapy with this hand-held control unit:

Use the intensity and stimulation frequency buttons to customise the effect of your NAIS THERAPY muscle and nerve stimulator. The most suitable **activity patterns** (i.e. therapy stimuli) and **stimulation frequencies** have been pre-programmed in the form of five fully automatic programs. A flashing light indicates which setting is currently in use.

NAIS has developed 5 fully automatic programs for the following areas of application:

1.  neck and/or shoulder
2.  lower back
3.  legs, including thighs and calves
4.  arms
5.  joints (suitable for both arms and legs)

ON/OFF + Intensity

This switch is for turning the unit on and for controlling the **stimulation intensity** (1-10). In other words, the therapy can be more or less intensive to suit your personal requirements. When the appliance has turned itself off automatically, reset the switch to "0". The unit works again when you turn it back on.

Stimulation Frequency

(only works in custom mode)
Nerve tissue can be stimulated electrically across a broad frequency spectrum. The stimulation frequency button enables you to set **the rhythm of the therapy stimuli**, on a scale from 1-10. This is important, as the response of nerve fibres can vary, depending on the sickness or injury, or any medication you may be taking.

The stimulation frequency cannot be customised when using a fully automatic program setting.

Impulse Modulation

i. e. the **continuous changing of activity patterns** (intensity of therapy stimulation).

You can add this effect by pressing the yellow button. The

The NAIS pads can also be controlled manually, i.e. without using fully automatic programs (*custom mode*).


Press

Kneading Impulse or **Throbbing Impulse** or **Vibration Impulse**
to obtain anything between a gently relaxing and a highly intensive massage effect.

The display

The display on the hand-held control unit tells you ...

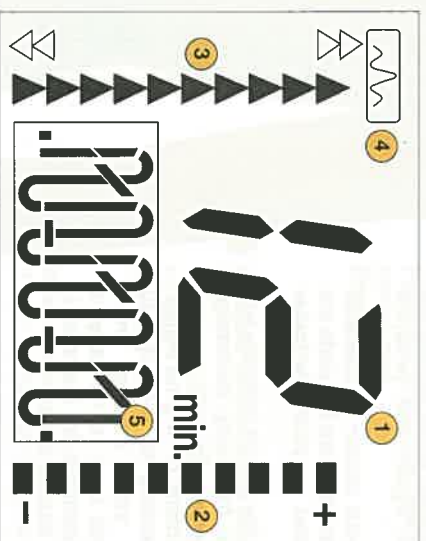How to position your NAIS pads - suggested applications:

Application notes:

For best effect, place the NAIS pads wherever the pain is at its worst. However, the painful zone may also be situated between two pads. To find out, feel carefully until you have located the points of pain.

You may need to adjust the positioning of the pads after a while. If the effect of the two pads differs (e.g. if therapy stimulation is less for one pad), change the position of the weaker pad. In this case make sure you first turn the hand-held control unit off and then restart it again when you have repositioned the pad.

The NAIS THERAPY muscle and nerve stimulator turns itself off as soon as you remove one of the NAIS pads from your skin (NAIS safety function). When you put the pad back again, the program restarts.

Consult your doctor or therapist when in doubt where the pads are best placed in your particular case.

Start with the lowest intensity and stimulation frequency, increasing the level gradually until you can feel a throbbing or tingling sensation. Depending on the selected intensity, your muscles may twitch. This level of intensity is used, for example, for the rehabilitation of muscle action.


Two programs are available for treating tense back muscles. Program no. 1 is for your upper section (neck and shoulders) and program no. 2 for your lower section (small of your back and lumbar region).

The illustration shows you a number of examples for back treatment:


- 1 Painful muscular tension in your back, caused through wrong sitting postures, etc.**
Place the NAIS pads on the pressure spots which, typically, have become hardened.
(Program no. 1)
- 2 Muscular tension at the top of your back / in your shoulders**
Place the NAIS pads directly on the painful zone.
(Program no. 1)
- 3 Muscular tension in the small of your back**
Place the NAIS pads on the most painful areas.
(Program no. 2)

- 4 Sciatic pain**
If your sciatic pain is due to muscular tension set program no. 2 Position


General tension in neck and shoulder muscles
(Program no. 1)


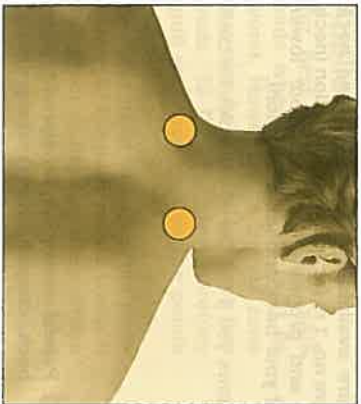
Muscular tension in your shoulders - caused, for example, by working at a computer
May also be placed on left shoulder, as an alternative.
(Program no. 1)

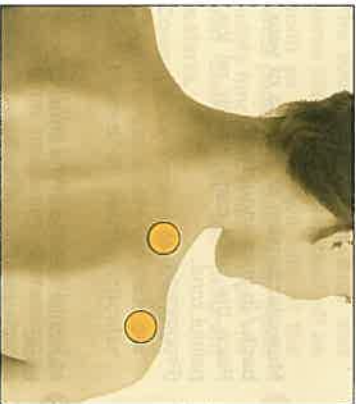Muscular tension in your neck and shoulders, caused by fatigue

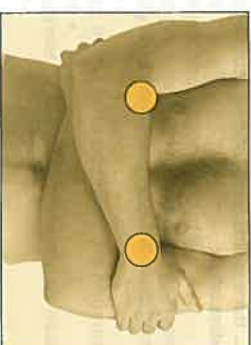Pain in the arm
The two illustrations show you ways of placing your NAIS pads on an arm. Muscular strains, circulation problems and/or the problematic healing of wounds in an arm can usually be treated with program no. 4.

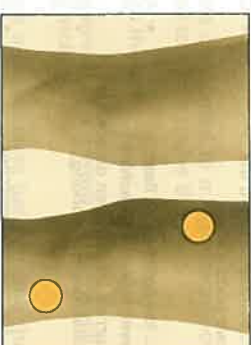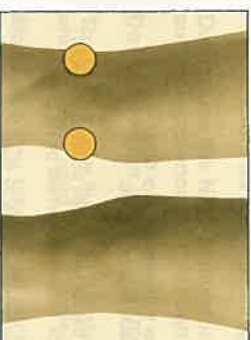Arm and knee joints
The illustrations show different ways of placing the pads on a knee. Program no. 5 is specially designed for stimulating painful joints.

Muscular strains, spasms and/or circulation problems

Legs can serve as an example for ways of treating muscular strains, etc. Place the NAIS pads lengthwise above and below the affected zone. Always use the special leg program no. 3 for legs. If you have varicose veins, consult your doctor before using NAIS THERAPY.

... there is little or no stimulation:

- Are the NAIS pads connected correctly to the hand-held control unit?
 - Check the push contacts between the pad and the cable.
 - Is the cable connected to the hand-held control unit?

- Have the NAIS pads been positioned correctly?
 - See pages 24-27.

- Are the batteries weak or flat?
 - Replace batteries.

- Have the batteries been inserted correctly?
 - Check correct placement of +/- poles.

... stimulation from the two pads differs:

- Change the position of the weaker pad.
- Are the batteries weak or flat?
 - Replace with new ones.

... the lights don't flash, even when the machine is turned off and then on again (machine not working).

.. the machine switches itself off automatically before the 15 minutes are over:

- Are the batteries weak or flat?
 - Replace with new ones.
- Is the treatment too long?
 - Shorten time of treatment.

... your skin gets irritated:

- Are the NAIS pads too dry?
 - Moisten with a drop of water. If this doesn't help or if the pads have lost their adhesion, replace with new ones.

- Have the pads been positioned correctly?
 - Place pads flat on the relevant skin section.

The NAIS muscle and nerve stimulator must not be used

- by patients with pace-makers and other implants,
- on babies and persons who are unable to express themselves,
- while carrying out an activity that demands your full attention, e.g. driving a car,
- while your body is wet, e.g. in the bath or under the shower.

Never position the NAIS pads in the areas of the skull or heart, or near the neck.

While in hospital, only use NAIS THERAPY after consulting your doctor.

We would advise against using NAIS pads during pregnancy, as no information has so far been gathered in this area.

If you have varicose veins, only use the NAIS stimulator on legs after consulting with your doctor.

Do not use the NAIS stimulator while connected to a high-frequency surgical appliance, as this might lead to burns under the pads.

Do not use after applying ointment to the skin.



The output may fluctuate if a short-wave or micro-wave therapy appliance is being used in the immediate vicinity.

Make sure you take off all necklaces and bracelets before using the NAIS stimulator.

The NAIS THERAPY stimulator is not a toy and should therefore not be accessible by children.

Never store or use under extreme conditions.

1. To open the battery compartment at the back of the hand-held control unit, shift the lid in the direction of the arrow, applying gentle pressure.
2. Insert alkaline batteries, 2 x 1.5 V (LR03/micro), making sure their poles face the right way.
3. Fresh, brand-quality batteries will last for approx. 200 x fifteen-minute treatments.
4. Rechargeable batteries are not suitable.
5. If the voltage of the batteries is too low, the appliance will turn itself off. Low voltage or output may mean that the NAIS stimulator suddenly stops.
6. If you do not use the NAIS stimulator for some time, you should take the batteries out of the hand-held control unit, to prevent leakage.

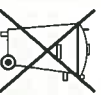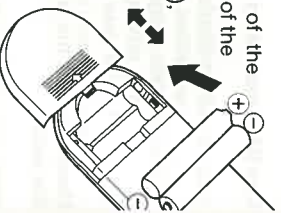Please help protect our environment by respecting national and/or local recycling regulations when disposing of the device and the batteries at the end of their useful lives.

Maintenance / cleaning the NAIS pads

The NAIS pads will last longer if your skin is clean and free from grease. If this is the case, they can be used approx. 20 times. When the pads lose their adhesion, you can extend their life to a total of 30 sessions by carefully dabbing them with a moist cloth. Always store the NAIS pads on the appropriate holder for the pads and the cable.

Spare pads can be obtained from specialist retailers, e.g. dispensing chemists and sanitary appliance shops. Make sure you only buy original NAIS accessories.

Infinitely adjustable Output voltage: 0 - 64 ± 10 V
0 - 64 ± 10 V
0 - 44 ± 10 V
200 µs

Length of impulse: 200 µs
Impulse polarity: automatic alternation

Infinitely adjustable stimulation frequency:

- Throbbing impulse 1 - 20 Hz
 - Kneading impulse 50 Hz
 - Vibration impulse 4 - 30 Hz
- length of burst 0.5 - 2.5 sec.
length of pause 0.5 - 2.5 sec.
length of burst 0.75 - 4.0 sec.
length of pause 0.5 - 1.0 sec.

The technical data refer to a 500 Ohm load resistor.

Programs: 15 minutes of alternating activity patterns (therapy stimuli) for each program

1. Neck/shoulders 4 cycles @ approx. 210 sec.
2. Lower back 5 cycles @ approx. 150 sec.
3. Legs 5 cycles @ approx. 150 sec.
4. Arms 5 cycles @ approx. 150 sec.
5. Joints 5 cycles @ approx. 150 sec.

Turns itself off automatically after 15 minutes.

Batteries: 2 x 1.5 alkaline (LR03/micro)

Life of batteries: approx. 200 treatments @ 15 minutes

Battery check: If the batteries are weak, the machine turns itself off, and they need to be replaced.

Classification: BF

Transport and storage conditions: Ambient temperature -10°C to +60°C
Relative humidity 10% to 95%
Operating conditions: Ambient temperature +10°C to +40°C
Relative humidity 20% to 85%

Example program no. 4 = arms	
Throbbing/slow 15 sec.	5 cycles
Kneading/fast 30 sec.	@ approx. 150 sec.
Kneading/slow 30 sec.	
Vibration/fast 15 sec.	
Vibration/slow 15 sec.	

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GB Warranty

This device has been manufactured with all due care and subjected to stringent testing before leaving our factory.

The NAIS THERAPY Muscle and Nerve Stimulator is guaranteed for 12 months as from the date of purchase as evidenced by the invoice and/or warranty card bearing the dealer's stamp.

During this warranty period, NAIS Wellnesslife GmbH will replace all defective parts free of charge, including the labour costs incurred for the warranted repairs. The warranty period is not extended by such repairs. Replaced parts become our property. All other rights, particularly the right of rescission, reduction in price or compensation, are excluded as a result of this warranty. The warranty does not cover defects due to misuse, non-compliance with the Instruction Manual or unauthorized modifications. Batteries and wear parts are not covered by this warranty.

If the complaint is not justified, we shall be entitled to demand an adequate fee, cash on delivery, for testing and mailing the device.

Please send the complete device with the stamped warranty card and/or purchase invoice to the NAIS customer service or the authorized after-sales service centre in your country (see overview) for any repair or warranty performance. In case there is no NAIS customer service centre in your country or territory or you are not aware of its address, please contact the outlet where you purchased the device.

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Cardiacare Ltd.
Wimed S.p.A.
Blokkandpack B.V.