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### ULTRAFAST THERMOMETER FOR BABIES

Before use, read the instructions and warnings given in this handbook carefully and keep the handbook in a safe place. Digital thermometer with a large LCD display for easy reading. Fast and accurate fitted with a sound signal when the temperature reading has been taken. Stores the last read temperature, to evaluate any variations. Extremely easy to use, it is kept in a practical case to avoid risk of breaking. The flexible probes can be washed, and it turns off automatically.

#### GENERAL WARNINGS

- 1 - Self-measuring means control, not diagnosis and treatment, any unusual readings must be discussed with your doctor. Under no circumstances must you change the medicine doses that have been prescribed by your doctor.
- 2 - Do not squeeze, bend, drop or break the thermometer.
- 3 - Keep it out of the reach of children.
- 4 - Do not expose to sunrays or direct sources of heat. Keep in a dry clean place.
- 5 - Do not immerse the thermometer in water or other liquids, only the probe is waterproof.
- 6 - This thermometer reads only the armpit, rectal or oral temperature. Do not use it for reading auricular temperatures.

#### INSERTING/REPLACING THE BATTERIES


If the display shows the battery signal, it means they need replacing. 1 - Pull the battery cover in the direction shown on the arrow. 2 - Use a pointed object to remove the old batteries. 3 - Place the new batteries with the positive pole (+) pointed upwards and the negative pole pointed downwards. 4 - Replace the cover. 5 - Never dispose of batteries with ordinary rubbish.

#### DESCRIPTION OF APPLIANCE

#### INSTRUCTIONS FOR USE

1 - Press the ON/OFF key. For a few seconds all the function symbols light up on the display. After the thermometer displays 36.5°C. When the °C symbol starts flashing, the thermometer is ready for use. 2 - Place the probe as follows: **IN THE MOUTH:** place the probe under the tongue and keep the mouth closed for about 10-15 seconds, until the temperature stabilises. The standard temperature for this type of reading is 37°C. **UNDER THE ARM:** place the probe under the armpit and keep the arm well pressed against the chest for about 10-15 seconds, until the temperature stabilises. The standard temperature for this type of reading is 36.7°C. **IN THE RECTUM:** lubricate the probe with a water-soluble gelatinous solution. Do not use Vaseline. Insert the probe into the rectum for about 1 cm., until the temperature stabilises (about 10-15 seconds). The standard temperature for this type of reading is 37.0°C.
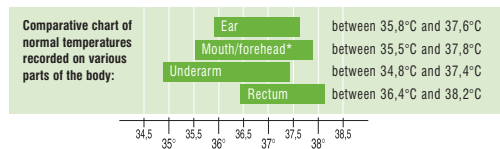
### ULTRAFAST THERMOMETER FOR BABIES

**3** - When the thermometer gives off the sound signal, remove it. The reading is terminated and shown on the display. The thermometer turns off automatically after 10 minutes without use (in the meantime the reading is still displayed) otherwise you can turn it off by pressing the ON/OFF button. If the temperature is above 44°C, the display shows HC, if it is below 32°C, the temperature shows LC. 0434 (Dir. 93/42/EEC)

#### MEMORY FUNCTION

The thermometer stores the last reading in its memory, to evaluate any variations. To recall the last reading, press the ON/OFF button for a few seconds until the reading appears on the display; in the bottom right corner of the display, "M" appears. Release the ON/OFF button to continue with measuring the actual temperature.

The thermometer stores the memory of last 5 reading, to evaluate any variation. To recall the last 5 reading, keep pressing the ON/OFF button and LCD displays 188.8 [°], following by indication 1M, 2M-5M while pressing the button each time.

#### CHANGING THE UNIT OF MEASUREMENT

Press the on/off button. The screen shows 188.8 ( ) and then 1M. When the screen displays ( ) and once the M symbol disappears, select the unit of measurement °C or °F by pressing the button. When the symbol °C or °F flashes, the thermometer is ready for use.

#### MAINTENANCE

Clean the thermometer probe with a cloth dampened with a 0.02% hypochlorite sodium solution or with alcohol or with any type of topical disinfectant before and/or after each use. The rest of the thermometer must be cleaned with a dry cloth. Keep the thermometer safe in its special case.

**Do not throw old batteries away. Place them in a container specially provided for worn batteries so that they may be collected and recycled.**

#### TECHNICAL DATA

Measuring range: 32°C-44°C  
 Tolerance: +/- 0.1°C 35°C-39.0°C, otherwise +/- 0.2°C  
 Accuracy: 0.1°C  
 Working conditions: 0°C to 50°C, Humidity: 10-83%  
 Conservation conditions: -20°C to 60°C  
 Humidity: 10-83%  
 Battery: 1.55Vdc manganese alkaline, SR41, LR41, UCC192  
 Battery autonomy: approx. 100 hours continuous use  
 Size: 126.5 mm x 18 mm x 10 mm  
 Weight: 10 g (with battery)

#### MANUFACTURER

**VALEO CORPORATION**  
 Formosa Office Building  
 4F-3, 161, Sung-The Road  
 Taipei, Taiwan, R.O.C.

#### WARRANTY

This product has a 2-year warranty against defective material and workmanship. During this period, such defects shall be repaired free of charge (proof of purchase should be presented whenever you have a complaint during the warranty period). This warranty does not cover damage caused by accidents, misuse or negligence. In case of complaint, please contact the shop from which you bought the appliance before contacting anyone else.

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