## INSERTION DES PILES

- INSERTION OF BATTERIES • ENLEGEN DER BATTERIEN • INSERIMENTO DELLE PILE**  
**• INSETION DE LAS PILAS • INSETTING AV BATTERIENE • INSERAREA BATERIILOR**  
**• INSEÇÃO DAS PILHAS • HET PLATSEN VAN DE BATTERIEN • WKŁADANIE BATERII**

**x 3**  
**AA LR6**  
**(1,5V)**

- non fournies • not supplied
- nicht mitgeliefert • non fornite
- no incluidas • følger ikke med
- nefurnizate • não fornecidas
- niet meegeleverd • poza zestawem

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

## Description of the unit BODY CONTROL DUO MINCEUR or BODY CONTROL DUO MINCEUR

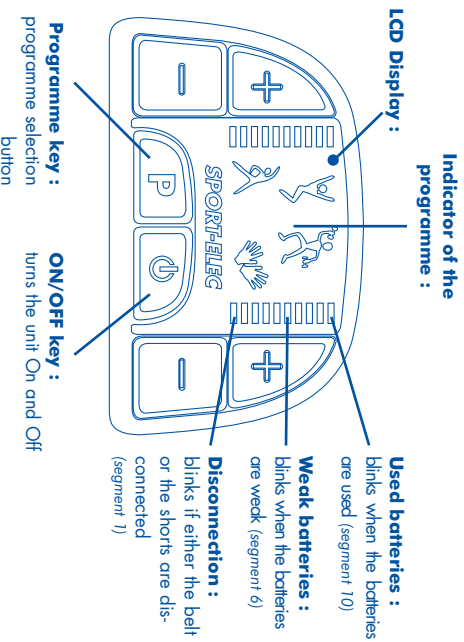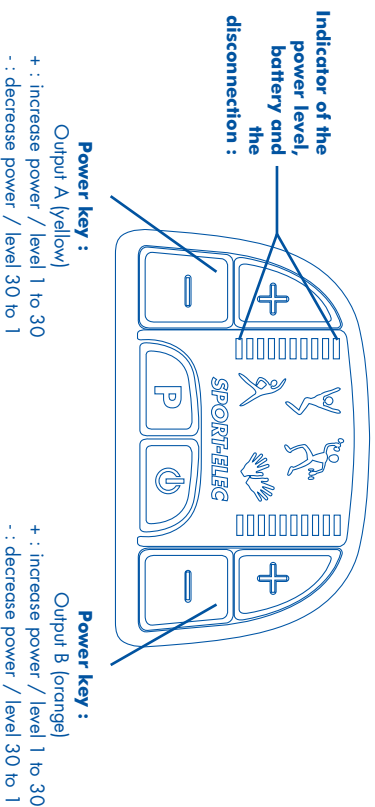NB : the segments 1, 6 et 10 have a double function.

- 11 -

## Programs

- PG 01 Firmness**  
20 minutes to warm up and to tone your muscles.  
Total duration : 20 min
- PG 02 Muscle care**  
40 minutes of body building to firm up, but also to increase muscular volume after a period of inactivity.  
Total duration : 40 min
- PG 03 Muscle Force**  
20 minutes of more intensive musculation to work the target areas.  
Total duration : 20 min
- PG 04 Massage**  
20 minutes to soothe tired muscles and allow you to feel the well being.  
Total duration : 20 min

### MUSCULAR KEEP FIT DEVICE SPORT-ELEC®

You have recently acquired the electro-musculation device with accessories, designed, produced and manufactured by **SPORT-ELEC Institut** (France). We thank you for your purchase and suggest that you read attentively the following paragraphs for an optimal use of the device.

**SPORT-ELEC** is specifically designed to enhance the abdominal area and the buttocks. The device has 79 frequency variations over 4 programmes. The 2 generators allow the optimal distribution of muscular contractions and make this device a highly effective product.

### WARNING

#### READ THE FOLLOWING INFORMATION BEFORE USE :

The generator is equipped with a safety system to avoid incorrect handling of the device by users.

Prior to use, check the following points :

- the device must be connected to an accessory placed on the body (if you turn on a device while no accessory is connected and press the "+- key to increase the power, the device will stop (safety stop).
- Apply the electro-conducting cream to the electrodes (supplied with the device).
- The device must be fully connected to all the 4 connecting plugs of the accessory.
- Adjust accessories carefully to ensure contact between electrodes and the skin.

**CAUTION:** If you move with accessories attached and a partial loss of contact between electrodes and the skin occurs, the device will stop (safety stop).

**CUSTOMER SUPPORT + 33 (0)2 32 96 50 50**

- 12 -

## I. INSERTION OF BATTERIES

Check that the apparatus is switched off. Open the "batteries" compartment, insert three 1,5 V AA LR6 batteries into the compartment then to close the case. Comply with + and - indications (**photo II - A page 1**). For an optimal use of your apparatus, it is advised to use good quality batteries.

## II. STARTING THE UNIT

### A / PRESENTATION OF THE UNIT

Your electrostimulation device is an appliance in the point of the technique, it respects the current standards in the field of medical application. The interface "man machine" is simplified and realized by means of pictograms and function keys. One of four, the pictogram in black indicates you the program chosen by the "P" key.

Two scales (1 by output) indicate you 30 power levels.

You also find 6 order keys (**scheme, page 3**) :

- 1 ON/OFF (red key ) - a short press
  - 1 "P" (grey key) to select the programs - a short press
  - 2 "+" / "-" yellow keys and 2 "+" / "-" orange keys to adjust the intensity of the stimulation (a short press).
- A long (1 sec.) press on the "-" key allows to return the power level rapidly to zero.

### B / B / UNIT FUNCTION

#### STEP 1 - Disconnection

Your device is endowed with an accessories disconnection detector.

If there is no contact with the human body and the accessory, the power is stopped on the concerned output (channel A = yellow / channel B = orange).

You are warned of a disconnection by blinking of a segment during 30 sec., after 30 sec. the blinking stops (**diagram, page 3**).

If no accessory is in touch with the body, the power cannot be emitted.

#### STEP 2 - Batteries

- If the batteries are "good", there is no blinking.
  - If the batteries are "weak", two segments (the 6th segment) blink (**scheme page 3**).
- Your device still works.
- If the batteries are "used", two segments (the 10th segment) blink (**scheme page 3**).
- Your device stops, even in the course of a program. Change the batteries.

#### STEP 3 - ON / OFF Switch "

- A short press on the "
- Every time under power, the device makes an auto-test on 1 essential parameters at the level of the safety of its operation. During the auto-test segments and pictograms are lit during 0,5 second.
- If an internal defect in the device is detected, all the pictograms flash. Send back the device to the After Sale Service.

- 13 -

- The auto-test does not take into consideration the used batteries. When the batteries are used the auto-test cannot be performed.
- One long press (over 3 sec) on the " key stops up the auto-test. In this case, stop the device by pressing the " key and then press shortly the " key again.

#### STEP 4 - Choice of a program "P"

It is enough to make one or several presses on the "P" key, the programs ravel. The programs are announced by pictograms. As soon as a program is selected, corresponding pictogram lights (**see page 3**).

Once the program is chosen, you can start it by increasing the power on the only of power keys "+". During the program the key "P" is inactive.

#### STEP 5 - End of program - Stop

- At the end of the session the device stops automatically.
- To stop the device in the course of program, you just have to turn 2 power keys to zero or to press the " key.
- In case of emergency, press the " key.

#### STEP 6 - Adjustment the power

You have 30 levels from 0 to 30 to adjust the power (there are three power levels between two segments).

## III. ELECTRO CONDUCTING CREAM

**SPORT-ELEC®** cream is used with accessories. It may be used for a manual massage of the body as well.

The cream facilitates the contact between electrodes and skin. Its massaging and electro-conducting properties improve comfort and efficiency of the stimulation.

If you wish to buy again some cream, order it directly at 02 32 96 50 50, on our website [www.sport-elec.com](http://www.sport-elec.com), or contact your local retailer.

- **Do not use the belt or the shorts without cream.**
- **Do not use any other cream except the one supplied by SPORT ELEC® Institut**

#### SAFETY PRECAUTIONS BY USING THE CREAM

- In case of allergy, stop to apply the cream.
- Do not swallow.
- Do not leave with the range children.
- In case of contact with eyes rinse with water.

## IV. BODY CONTROL MINCEUR

### A / BELT FUNCTION

- 14 -

**PRESENTATION OF THE BELT :**

- **Front side : (photo I - F page 2)**
- 1 - Wire storage
- 2 - Connection plugs
- 3 - Scratch strap
- **Electrode side : (photo I - D page 2)**
- 4 built-in electrodes

**STEP 1 : Connection and installation of the belt.**

- 1 - Make sure that your **SPORT-ELEC®** device is switched OFF.
- 2 - Take the connectors of the unit and plug them to the belt plugs by 2 groups **(photo I - F page 2)**.
- 3 - Apply a thick layer of the electro-conducting massage cream on your abdominal area and the on the built-in electrodes of the belt **(the convex part - photo I - E page 2)**
- 4 - Place the belt on your abdominal muscles in such a way that the electrodes should be centred on your abs. Adjust them using the "stretch" straps making sure that the built-in electrodes are in contact with the skin **(photos I page 2 : G - abdominal muscles)**.

**STEP 2 : Session start**

Once your device connected, make yourself comfortably and start your device following the procedure by "**Starting the unit**" detailed previously in paragraph II.

Do not try to increase the intensity too quickly. If the power seems to you too strong, decrease it by pressing briefly on the "r" key. A long press on the "r" key returns the power to its minimal level. This operation does not switch off the device.

Do not hesitate to place the belt slightly to the left or right (after switching the power key on zero) in order to feel the contractions on the whole area of abdominals.

**If an electrode hurts or if you feel some tingling, it is necessary to stop the device and to verify whether the electrode is well covered with cream.**

**STEP 3 : End of the session**

Once the programme has ended, switch off the device and unplug it from the belt. Remove the surplus of cream from the built-in electrodes with a moistened cloth and dry the belt. Store carefully your material.

## V. BODY CONTROL DUO MINGEUR

### A / BELT FUNCTION :

See the previous paragraph "BELT FUNCTION".

### B / SHORTS FUNCTION :

**PRESENTATION OF THE SHORTS :**

- **Front side : (photo I - F page 2)**
- 1 - Wire storage
- 2 - Connection plugs
- 3 - Scratch strap
- **Electrode side : (photo I - D page 2)**
- 4 built-in electrodes

**STEP 1 : Connection and installation of the shorts.**

- 1 - Make sure that your **SPORT-ELEC®** device is switched OFF.
- 2 - Take the connectors of the unit and plug them to the shorts plugs by 2 groups **(photo I - F page 2)**.
- 3 - Apply a thick layer of the electro-conducting massage cream on your buttocks area and the on the built-in electrodes of the shorts **(the convex part - photo I - E page 2)**
- 4 - Place the shorts on your buttocks **(photos I page 2 : G - buttocks)**.

**STEP 2 : Session start**

Once your device connected, make yourself comfortably and start your device following the procedure by "**Starting the unit**" detailed previously in paragraph II.

Do not try to increase the intensity too quickly. If the power seems to you too strong, decrease it by pressing briefly on the "r" key. A long press on the "r" key returns the power to its minimal level. This operation does not switch off the device.

Do not hesitate to place the shorts slightly to the left or right (after switching the power key on zero) in order to feel the contractions on the whole area of buttocks.

**If an electrode hurts or if you feel some tingling, it is necessary to stop the device and to verify whether the electrode is well covered with cream.**

**STEP 3 : End of the session**

Once the programme has ended, switch off the device and unplug it from the shorts. Remove the surplus of cream from the built-in electrodes with a moistened cloth and dry the belt. Store carefully your material.

**CONTRA-INDICATIONS**  
**Do not use SPORT-ELEC® :**

- \* on the heart area
- \* If you carry a pacemaker
- \* If you are pregnant
- \* during first 6 weeks after baby birth
- \* In case of neuromuscular diseases
- \* In case of hemorrhagic diseases

## VI. Recommendation of use

The unit is recommended for physiotherapists for functional re-education.

### WARNINGS :

Do not use the device by driving or by working on a machine.  
Do not use this electro stimulation unit in a damp environment, a bathroom, a sauna, or close to electric machinery connected to mains and earthed or even piping earthed.  
The simultaneous connection of a patient to a high frequency surgical unit may cause burns at the points of contact of the stimulator's electrodes, and the stimulator may be damaged. Operation in close proximity of a short wave therapy unit may cause instability in the stimulator's output power.

### RECOMMENDATIONS :

Please apply regulations enforced in your country when disposing your product at its end of life. Therefore, please discard it in a location planned to this effect, ensuring it is disposed in an environment-friendly manner

### PRECAUTIONS OF USE

- \* Keep the belt and the shorts away from heat
- \* Do not leave within reach of children
- \* **Ask your doctor's advice :**
- if you have electronic or metallic implants (pins, etc.)
- in case of skin condition (wounds,...)

- Use only on those parts of the body for which the accessories were designed.
- Functions only on batteries
- Not to connect to another source of power

### RECOMMENDATIONS OF USE

- You can restart a program's cycle on another muscle zone or change the program while remaining on the same muscle zone, without exceeding however 40 min use per muscle zone and per day.
- For optimum use, refer to the diagrams shown in this leaflet.
- Clean the unit and its cables with a cloth dampened with a disinfecting solution (mercury type) diluted in water at 1/10th.

### TIPS

#### EXPLANATIONS FOR SMALL SIZES :

If the belt overlaps with electrodes on the body, hold the overlapping part of the belt by means of scratch sections.

### CLEANING AND MAINTENANCE CONDITIONS FOR THE BELT AND SHORTS.

- CAUTION : Machine washing and soaking is PROHIBITED.
- Do not dry-clean
- Do not use bleach
- Clean accessories with a humid sponge soaked with detergent and wipe with a cloth.
- Clean the unit with a cloth soaked in with disinfectant (Mercury type) diluted with water at 1/10th.
- Clean the electrodes with a mix of water/hypochlorite detergent, then rinse thoroughly.

### SAFETY CHECK

Every 24 months, check the following points :

- 1 - Check for superficial degradation of the equipment
- 2 - Check proper operation of the unit :
  - \* legible display
  - \* switch, contact
- 3 - Measure electric values :
  - \* with a 500-ohm load, the current must not exceed 100 V.
- 4 - Check reliability of accessories :
  - \* adhesive electrodes in good condition (no tears, no material detests)
  - \* wire of adhesive electrodes in good condition.

## VII. TECHNICAL SPECIFICATIONS

UNIT :	
Type of current	Dual phase
Number of programmes	4 programmes
Number of phases variations	79
Number of adjustable independent outputs	2
Current	From 0 to 70 mA - From 0 to 1000 Ω
Frequency range	1 to 100 Hz
Pulse width range	160 to 320 µs
Power	3 x 1,5 V batteries type LR6 (not provided)
Unit dimensions and weight	69 x 43 x 87 (h) mm - 106 g
Automatic stop at end of each program	Yes - about 4 min 30 sec
Activation safety	100%
Detection absence of electrodes or of contact	yes detection independent on each of the 2 outputs
Indicator of power	bar graph - LCD display
Indicator of weak and/or defective batterie	yes - LCD display
Indicator absence contact belt / shorts	yes - LCD display
Indicator of storage	picograms - LCD display
Relative humidity of storage	0°C to + 45°C
Temperature of use	10% to 90%
Relative humidity of use	+ 5°C to + 45°C
Connection device belt / shorts	20% to 65%
Caution : Read the manuel thoroughly	wires - plugs ø 2 mm mole
RoHS : The unit are free of materials likely to have a negative impact on the environment	<b>RoHS</b> (2002/95/EC)
Guarantee	2 years
<b>STANDARDS</b>	<b>CE</b> (Directive 89/336/CEE)
• Unit : EN60601-1 EN60601-2-10	
• Belt / shorts: ISO 10993	
• Electro-conducting massage cream 150 ml / 5,1 Fl Oz : 76/768/CEE	
• Classification type BF	
<b>BELT AND SHORTS :</b>	<b>GUARANTEE :</b>
• Size (unisex) : from 34 to 60	The belt, the shorts and the cream are designed for personal use. For reasons of hygiene, there will be neither trade-in nor exchange.
• Length of the unfolded belt : 126 cm	The belt life is limited to 10 months.
• Fabric : polyester	
• Built-in electrodes : carbonated silicon.	

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