**EKS**

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**GB** The new generation of EKS Body Analyzer Scales calculates your body composition using unique **TAMM "Triple Aided Multi-frequency Measurement" technology**. This technology is superior to basic bioelectrical impedance analysis (BIA) technology and will provide you with highly accurate readings.

With the basic BIA, the individual resistance in the body tissue is measured by a low-level electrical signal that passes through your body. It is easy for the signal to flow through fluids in the muscle and other tissues but it meets resistance as it passes through body fat.

Different frequency signals determine the proportion of body fat, water, muscles and bone mass in the human body. As the sizes and the distance between body cells vary, they react differently to different frequencies.

Contrary to the usually applied 50 kHz signal to determine body fat and water percentage, EKS clinically validated studies show that 54 kHz is the optimum frequency for a more precise measurement.

In order to guarantee the precision of the measurement of all relevant tissues with different characteristics, two further frequencies 52.6 kHz and 56 kHz run parallel with and support the 54 kHz. This is the meaning of the phrase "Triple Aided Multi-frequency Measurement".

EKS Body Analyzer Scales are designed for healthy living. Their unique technology enables the accurate measurement of key body health indicators and come equipped with an "athlete" mode for those involved in intense physical activity.

#### Disclaimer:

Warning! The following people should not use this analyzer:

- People with battery driven implants
- Pregnant women

This device is for domestic use only and is not suitable for children under 18.

For people trying to loose or gain weight, it is recommended that they do it in consultation with a doctor.

**EKS cannot be held responsible for any damage or injuries caused by an analyzer to the user or any third party.**

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If the scale has been programmed for the user, it will automatically recognize the user and show the correct user number. If the scale has not been programmed see Section 4 - Programming of personal data

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If this weight or level is programmed for the user, it will be recognized automatically and the correct user number will be shown. If this weight or level is not programmed, see section 4 - Programming of personal data.

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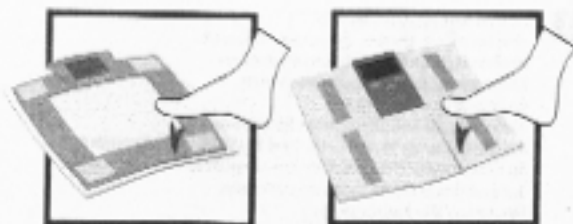
2-4 sec

2

Wait for "0.0KG" or "0.0LB" or "0:05T" in display before standing on scale.

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**GUARANTEE**

Before returning the scale, check that you have carefully followed all the instructions above. The guarantee period has been increased to 5 years (as from the date of purchase). This guarantee includes the repair fee or charge of any fault and the free replacement of any parts judged by us to be defective. The guarantee does not cover any damage resulting from the incorrect use of the scale. In the advantage of the guarantee, take the defective scale to the point of sale, together with the guarantee certificate bearing the retailer's stamp and signature. We reserve the right to make technical modifications.

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**CE CONFORMITY**

This scale satisfies all European requirements on spurious radiation and resistance to parasitic interference (in compliance with EEC directive 89/336). If interference occurs during the weighing, due to a strong source of electromagnetic radiation (cordless telephones, transmitting devices etc.) move the scale away from the source of emission or switch off the appliance responsible for the interference and repeat the weighing operation.

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**REGULATION FOR WASTE DISPOSAL 2002/96/CE**

This product is encompassed by the category for electric and electronic components, which covers products that may contain dangerous substances for human health and the environment. The crossed-out bin symbol means that the product is subject to strict disposal arrangements under European legislation (Directive 2002/96/CE). The user must dispose of the product either by taking it to a shop that sells similar products ("Obligation of return") or by putting it into the hands of a local refuse collection unit for selective recycling. Any product collected in this way will be dismantled, cleaned and disposed of in an environmentally safe way. Selective recycling helps to reduce the amount of waste and deals safely with any toxic substances. This decreases harmful effects on the environment in accordance with the Waste of Electric and Electronic Equipment (WEEE) Directive.

To obtain more information about waste disposal procedures, please contact your local refuse collection unit.

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