solac

MOD. ME7705

**MASAJEADOR SHIATSU
SHIATSU MASSAGER
MASSEUR SHIATSU
SHIATSU MASSAGEGERÄT
MASSAJADOR SHIATSU
MASSAGGIATORE SHIATSU
MASSAGE APPARAAT SHIATSU
SHIATSU MASÁŽNÍ PŘÍSTROJ
APARAT DO MASAŻU SHIATSU
SHIATSU MASÁŽNY PŘÍSTROJ
SHIATSU MASSÁZIROZÓGÉP
УРЕД ЗА ШИАЦУ МАСАЖ
SHIATSU MASAŽER
APARAT DE MASAJ SHIATSU**

**INSTRUCCIONES DE USO
INSTRUCTIONS FOR USE
MODE D'EMPLOI
GEBRAUCHSANLEITUNG
INSTRUÇÕES DE USO
ISTRUZIONI PER L'USO
GEBRUIKSAANWIJZING
NÁVOD K POUŽITÍ
INSTRUCJA OBSŁUGI
NÁVOD NA POUŽITIU
HASZNÁLATI UTASÍTÁS
инструкция за употреба
UPUTE ZA UPOTREBU
INSTRUCȚIUNI DE UTILIZARE**

5- 7 • ESPAÑOL
8-10 • ENGLISH
11-13 • FRANÇAIS
14-16 • DEUTSCH
17-19 • PORTUGUÊS
20-22 • ITALIANO
23-25 • NEDERLANDS
26-28 • ČESKY
29-31 • POLSKA
32-34 • SLOVENSKY
35-38 • MAGYAR
39-41 • българск
42-44 • HRVATSKA
45-47 • ROMANA

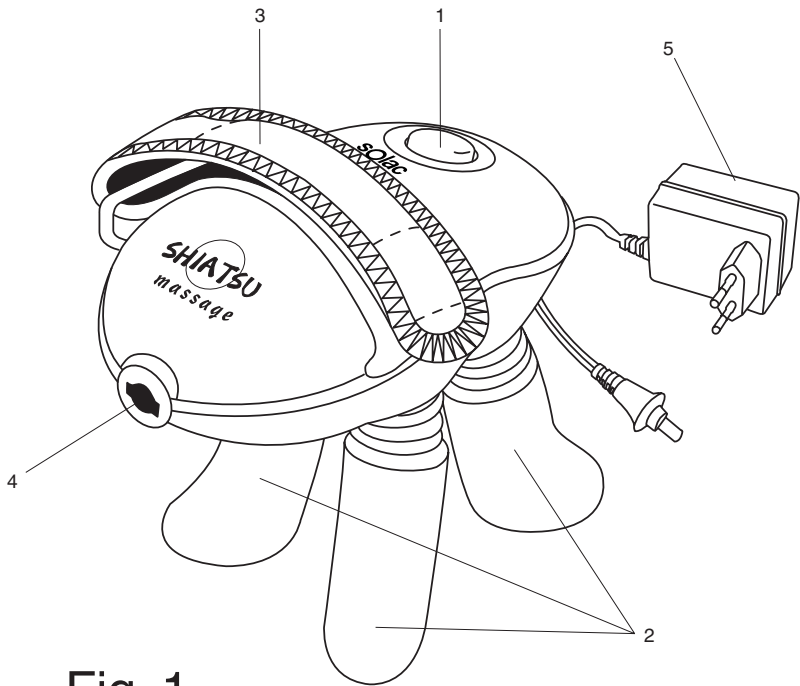
Fig. 1

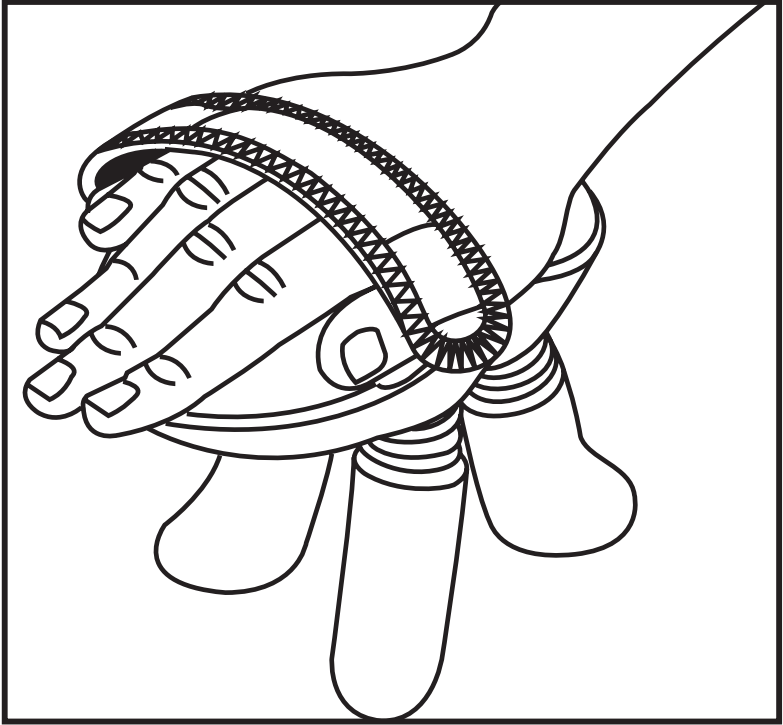
Fig. 2

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INTRODUCTION:

- The Shiatsu technique refers to the use of the fingers and palms of the hand to apply pressure to certain areas of the body to maintain and increase health. It is also a method for helping to cure specific diseases.
- Shiatsu stimulates the parasympathetic nervous system, which relaxes the body and activates the immunological system by mobilising defences and producing a feeling of wellbeing.
- The main therapeutic effects of this type of massage are as follows:
 - * Releases muscle tension
 - * Improves blood circulation
 - * Stimulates the circulation of body liquids
 - * Stimulates self-regulation of the body
 - * Treats various disorders: insomnia, stress, agitation, nervousness, headaches, migraine, neck pain, back pain, waist pain, joint pains, muscle tension...

IMPORTANT:

- Please read these instructions before switching on the appliance.
- Check that the voltage in your home matches that stated on the specifications plate on the appliance.
- Unplug the appliance when it is not in use. To unplug it, pull on the plug, not the cable.
- Do not use the appliance near a bath, shower or swimming pool. Do not use it while taking a bath or shower.
- Do not try to pick up the appliance if you accidentally drop it in water; unplug it immediately.
- Do not plug in or unplug the appliance while it is still in use.
- Keep it out of the reach of children and the disabled.
- Do not use in areas with cuts, burns, inflammation, rashes, ulcers, varicose veins or other disorders.
- Pregnant women and people with pacemakers, thrombosis or unexplained pain should consult their doctor before using the appliance.
- Do not use the machine in the abdominal area if you are in pain or on the carotid artery on the right or left of the neck.
- The massage should be pleasant and agreeable. If you notice pain or discomfort at any time, stop using the machine and consult your doctor.
- Do not use this appliance for purposes not given in these instructions.

- Unplug the appliance from the mains if it is to be unused for a long time
- To keep the appliance in a good state of repair, turn it off for 10 minutes after 20 minutes of uninterrupted use.
- If the mains cable is damaged, take the appliance to an Authorized Service Centre to have it replaced.

IMPORTANT: when you want to dispose of the appliance, **DO NOT** throw it in the rubbish **UNDER ANY CIRCUMSTANCES**. Take it to your nearest **CLEAN POINT** or waste collection centre. In this way, you will be helping to look after the environment.

MAIN COMPONENTS (FIG. 1):

1. ON/OFF switch (two speeds)
2. Flexible massage fingers
3. Adjustable hand strap
4. Adapter connection point
5. Adapter

CHARGING BEFORE USE:

- Turn the switch (1) off and plug the adapter (5) into the mains and to the connection point (4).
- When charging for the first time, leave the appliance on charge for 24 hours. After the initial charge, 11 hours will be enough to recharge the appliance. After a full charge, the appliance will work for 140 minutes.

MAINS CONNECTION:

- When the appliance has no charge, it can be used from the mains.
- Connect the adapter (5) to the connection point (4) and plug it into the mains.

OPERATION:

- This massager has four micro-motors which activate the flexible massage fingers (2) to imitate the movements of a Shiatsu masseur.
- Insert the hand between the strap (3) and the appliance, placing the palm of the hand on the appliance as shown in Fig. 2. Adjust the strap for maximum comfort.

- To turn the massager on, press the switch (1) once for minimum speed and again for a faster massage. To turn it off, press it a third time.
- When using the appliance, move the fingers on the area of the body to be massaged, holding the massager firmly.
- Do not massage the same area for more than 5 minutes. Massage each area for a short period of time and then move on to a different area.
- Use the massager for a maximum of 20 minutes for the appliance to last as long as possible.
- N.B. If you notice pain or discomfort, discontinue the massage immediately.

APPLICATION:

- We recommend you use the massager as follows, depending on the part of the body being treated:
- Back: the massage can be from top to bottom along the spine and in a circular motion in areas with greater muscle tension.
- Neck: apply the Shiatsu from top to bottom, from the base of the skull to the sides of the spine.
- Arms and legs: apply the Shiatsu with gentle taps from top to bottom along the length of the muscles.
- General massage: use the massager with circular movements all over the body.

MAINTENANCE AND CLEANING:

- The adapter has been designed exclusively for this product. Do not use any other kind of adapter.
- Always unplug the appliance before cleaning.
- Do not use abrasive cleaning products or liquids to clean the appliance. Use a damp cloth only.

BATTERY INFORMATION:

- The batteries should be removed only when they no longer work.
- All repairs and battery replacements should be carried out by an Authorized Service Centre.
- Do not throw used batteries away with your normal rubbish. Put them in the appropriate bins.

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