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© Irox Development Technology www.irox.com

Locus Electronics Limited
 17 Sandford Road
 Bexley heath Kent DA74AX U.K
 TEL: +44 7966 313 677

Cotronic Manufacturing
 Fl.4,Blk 4, Fu Yuan Ind Zone
 Jiu Wei, Xi Xiang, Bao An Zone
 People's Republic of China
 TEL : +852 2602 2021


CT120 Electronic Digital Thermometer adopts scientific technology of microcomputer. It features fast measurement, safety and simple operation.

Features

- Electronic digital thermometer
- Display of human body temperature
- Memory record (Last reading)


Keys & Display

1. LCD Display
2. Switch
3. Sensor

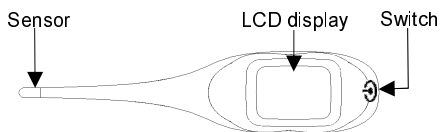Battery Installation

- 1) Open the battery cover with the use of a coin. Insert the coin at the back of the battery cover. Turn it counterclockwise to loosen it and take out the old battery.
- 2) Insert a new LR1130 battery. Put back the battery cover. Close it by turning it clockwise. Avoid damage to the waterproof ring.

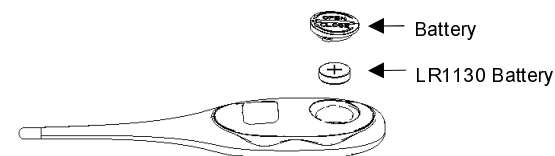How to use

- 1) Press the power switch. The thermometer starts a self-test routine and the screen shows the full segment with a sound of "beep", which indicates that the unit has sufficient power.
- 2) Memory function: After the full segment is shown, the last value with the letter "M" is displayed. Thereafter, the display shows Lo °C / Lo °F and the measurement begins.
- 3) During the measurement the "C" or "F" symbol is flashing and the temperature is rising.
Indications during the measurement:
a) "Lo °C" or "Lo °F": the temperature measured is below 32.0 °C / 90 °F.
b) Indication of the temperature in °C or °F.
c) "Hi °C" or "Lo °F": the temperature is above 43.0 °C / 109.9 °F.
- 4) When reaching the temperature maximum, a "beep" ton will sound for 10 seconds and the symbol "C" or "F" will stop flashing. The measurement will stop when the temperature did not rise for 16 seconds.
- 5) Indications after the measurement:
a) Normal temperature: "Beep" ton in a slow rhythm, the temperature is normal and is under 37.8 °C / 100 °F.
b) Fever alarm: "Beep" ton in a fast rhythm, the temperature is above 37.8 °C / 100 °F.
c) Temperature too high: "Beep" ton in a fast rhythm and the indication of "Hi °C" or "Hi °F", the temperature is above 43.0 °C / 109.9 °F.
d) Temperature too low: "Beep" ton in a slow rhythm and the indication of "Lo °C" or "Lo °F", the temperature is below 37.8 °C / 100 °F.
- 6) Turn off the unit by just slightly pressing the power switch. If the unit is left for about 9 minutes, it will automatically turn-off to save power.
- 7) When it is used to measure a child's temperature, close supervision should be provided by an adult.

Measurement:

You are advised to measure your temperature in morning.

- 1) Used in mouth: Please close your mouth for stability of temperature, breath calmly through nose. Deep breathing would affect the accuracy of measurement.
Normal temp.: About 37.0 °C / 98.6 °F
- 2) Used in anus or rectum: This is the most steady way to get the most accurate temperature. It is suited for both infants and children. Insert the sensor for about 2~3 cm. Vaseline can be used for easier insertion.
Normal temp.: About 37.6 °C / 99.7 °F
- 3) Used under armpit: Press and fix the thermometer with the arm. The temperature value of this measurement is not so accurate and not recommended for accurate measurement.
Normal temp. : About 36.7 °C / 98.1 °F
Measurement time: 3~5 minutes

Notes: The normal temperatures may differ from the different human bodies and measuring time.

Sterilization is required before every measurement.

Measurement is inaccurate when the battery is low-powered

Cleaning and maintenance:

The thermometer is waterproof and can be cleaned by water or sterilization lotion with soft cloth. Sterilization lotion below is recommended:

Name	Density
Medical alcohol	95%

Notes: Other lotion, density or cleaning method may damage the thermometer. The product doesn't need any calibration and is recommended to be used up to 3 years.

Remarks:

- 1) This CT120 Electronic Digital Thermometer is suitable for measuring human body temperature. It is a necessity in every household.
- 2) It is designed and manufactured in compliance with high standard requirements in accordance with EN12470-3. This thermometer performs with maximum accuracy and efficiency.
- 3) It has passed EMC emission test and immunity test in accordance with EN60601-1-2. It has passed biocompatibility test as well in accordance with MDD annex.
- 4)  Under the degree of protection against electric shock, it is classified to BF, in accordance with EN60601-1.
- 5) Be careful in handling it. Do not drop.
- 6) Never put it in boiling water for sterilization.
- 7) Do not put it close to any strong chemical solution and keep it away from straight sunlight or heat source. Avoid putting it near to strong electromagnetic source as this may cause interference in its performance.
- 8) It should be used in an environment with temperature of +10 °C ~ +35 °F (+50 °F ~ +95 °F) and with humidity of 40~80%. During transport, it must be stored in an environment with temperature of -10 °C ~ 40 °C (+14 °F ~ +104 °F) and with humidity of 30~90%.
- 9) Every 12 months or when the display becomes indistinct. When the symbol "▼" appears on LCD, the battery must be replaced to achieve accurate measurement.
- 10) Do not remove any parts except when changing the battery.
- 11) The scrap of product and battery are to be disposed according to national rules and regulations.
- 12) **Caution! Be sure the thermometer is kept away from children's reach.**

Specifications

Measuring Range:	32.0 °C ~ 42.9 °C (90 °F ~ 109.9 °F)
Measuring time:	About 3~5 minute
Accuracy:	±0.1 °C (±0.2 °F) for 35.0~39.0 °C (95 °F ~ 102 °F) & ±0.2 °C (±0.4 °F) for remaining range at an ambient of 22 °C (71.6 °F)
Display:	41/2 LCD digits (unit of display : 0.01 °C / 0.01 °F)
Battery:	LR1130 alkaline
Battery life:	About 250 hours or 1000 times.
Weight:	22 g

© Irox Development Technology www.irox.com

Locus Electronics Limited
17 Sandford Road
Bexley heath Kent DA74AX U.K
TEL: +44 7966 313 677

Cotronic Manufacturing
Fl.4, Bk 4, Fu Yuan Ind Zone
Jiu Wei, Xi Xiang, Bao An Zone
People's Republic of China
TEL : +852 2602 2021

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