Bilancia pesapersona
Bathroom Scale
Pese-personne
Personenwaage
Báscula

Biocheck

ITALIANO

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ENGLISH

- A) Battery compartment
- B) ON/OFF button for Body Fat Measurement (BFM) function
- C) —△ & ▽ selector for modifying height, age and sex settings in BFM mode
- D) MODE selector for switching among height, age, and sex in BFM mode
- E) ENTER button for confirming data in BFM mode
- F) RECALL button for recalling information stored in memory
- G) Unit of measurement selection button

USE OF THE APPLIANCE

This appliance must be used only for the use for which it is intended: as a bathroom scale and body fat measurement device. Any other use is incorrect and therefore hazardous. Polti S.p.A. shall not be held liable for any damage resulting from incorrect or erroneous use.

WARNINGS

Read the instructions carefully before using the appliance and keep them for future reference.

This appliance is intended exclusively for domestic use, and must not be used in hospitals or healthcare facilities.

Individuals with medical implants (cardiac electro-stimulators, etc.) should not use this device. If you have any doubts, ask your doctor for information.

Pregnant women should not use this device.

If the device is faulty or malfunctioning, do not attempt to dismantle it yourself, but contact an authorised Polti service centre.

If the appliance is dropped, have it inspected at an authorised service centre.

The appliance has no parts requiring maintenance other than the battery.

Wipe with a damp cloth after use. Do not use chemical products or submerge it in water.

Do not subject the appliance to excessive vibration or blows. Do not drop it onto the floor.

Do not leave the appliance where it will be exposed to the weather.

Do not leave the appliance near a heat source.

The appliance can measure using the metric system (kg/m) or the Imperial system (lb/in). Select the preferred scale using the sliding switch on the bottom of the appliance.



Always take used batteries to collection containers.

PRINCIPAL FEATURES

Your new BIOCHECK bathroom scale is designed to the highest quality standards.

It uses Bioelectric Impedance Analysis to determine the amount of body fat.

- High precision weighing technology
- Weighs up to 140 kg (308 Lb)

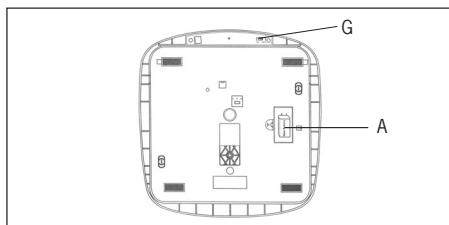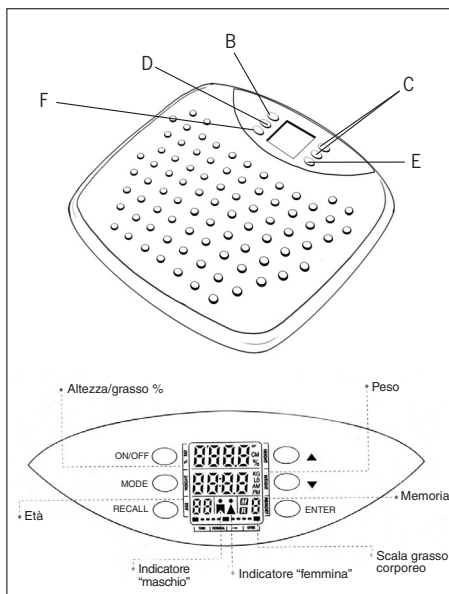- Scans every 100 g (0.2 Lb)
- Body fat scanning rate of 0.1% (within the 5% to 50% range)
- Height range 100.0 – 199.9 cm (3' 3.4" – 6' 6.7")
- Age range 10 – 80 years
- Measures both weight and body fat in about 10 seconds
- Can store data on ten people in memory
- Weight-only function comes on automatically when you get onto scale
- User may select kilograms/metres or pounds/inches
- Low battery indicator
- Runs on a 9v alkaline battery (6LR61 or equivalent), not included.

PREPARATION

Place a 9V alkaline battery (6LR61 or equivalent, not included) in the battery compartment (A), noting the + and – marks on the appliance and making sure it stays in the correct position. Select the desired unit of measurement (kilograms or pounds) using the sliding switch on the bottom of the appliance. Place the scale on a flat, stable surface to ensure accurate weighing.

USE

Set the scale the first time you use it. Gently press the scale's weighing surface. The display will show  and then the scale will automatically switch off. It is now ready for use.

Weighing

1. The scale comes on automatically; simply step onto the scale and wait briefly until your weight is displayed.
2. Your weight will appear on the liquid crystal display.
3. The scale will go off automatically after use.
4. Overload indicator. If the message "Err2" appears, IMMEDIATELY remove the weight from the scale or permanent damage may result.

Measuring body fat

Fig. 1

1. Press the (ON/OFF) button to start BFM mode; the previous setting will appear. The height indicator will flash (Fig.1)

Fig. 2

2. Press the (Δ) or (▽) keys to increase or decrease height.
3. Press (ENTER) to confirm the height setting. The age indicator will now flash.

Fig. 3

4. Press (Δ) or (▽) to modify age.
5. Press (ENTER) to confirm the age entered. The sex indicator will now flash.

Fig. 4

6. Press (Δ) or (▽) to modify the sex setting using the corresponding symbol.
7. Press (ENTER) to confirm the values you have entered. You will hear a long "beep".

Fig. 5

8. The  symbol will appear (Fig. 2)

Fig. 6

9. Step onto the scale with bare feet. Stand still while you are weighed. Your weight will be displayed for a few seconds (Fig. 3)

Fig. 7

10. The scale will now start measuring body fat. While this procedure is underway the message "BFM measurement in progress" will appear on the display (Fig. 4).

Fig. 8

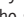
11. When BFM measurement has been completed, a long "beep" will sound. The appliance will display fat percentage, weight, age and sex for a few seconds and then go off automatically (Fig. 5).

Body fat measurement using the memory function

The appliance can store data on ten different people in memory. It can store height, weight, age and sex data for frequent use.

To store data in memory:

Repeat the procedure described in the section on "MEASURING BODY FAT" up to point 8. The memory symbol  will appear (Fig. 6). Press (Δ) & (▽) to select the memory position (0 – 9); the memory symbol will flash. If the selected position in memory is full, the symbol  will appear under the memory symbol (Fig. 7). Press (ENTER) to store the body fat measurement in the desired position in memory. A long "beep" will sound while memorisation is underway.

To recall information in memory

Press (ON/OFF) to start BFM mode. Then press (RECALL) to start up the memory recall function. The current setting will appear. Press (Δ) & (▽) to select the desired memory position (0 – 9). If the memory position is empty, no previous input will be displayed (Fig. 8). If the memory position is occupied, the information previously stored in memory will be displayed (Fig. 7). Press (ENTER) to accept the BFM settings stored in memory; a long "beep" will sound. Step onto the scale with bare feet. Stand still while you are weighed. Your weight will be displayed for a few seconds (Fig. 3). The scale will now start measuring body fat. While this procedure is underway the message "BFM measurement in progress" will appear on the display (Fig. 4). When BFM measurement has been completed, a long "beep" will sound. The appliance will display fat percentage, weight, age and sex for a few seconds and then go off automatically (Fig. 5).

Error messages

- Display Err0** Initialisation error. Try repeating the measurement.
- Display Err1** Incorrect use or body movement. Repeat the measurement.
- Display Err2** Overload indication. Remove the weight from the scale immediately or permanent damage may result.
- Display Err3** Body impedance is beyond technical limits (300 ohm – 1 kohm). Wipe the soles of your feet with a damp cloth and repeat the measurement.
- Display Err4** Body fat percentage falls outside the range (5% - 50%). Wipe the soles of your feet with a damp cloth and repeat the measurement.

Additional information

1. Ideal body fat percentages as recommended by experts:

Ideal body fat percentages	
Women	20% to 29%
Men	10% to 18%

ENGLISH

2. Consult your doctor to determine your ideal body fat percentage.
3. Always measure body fat at the same time of day (for instance, between 6 and 7 in the morning).

MAINTENANCE

Replace the battery when the display indicates "LO".

Remove the battery if you will not be using the scale for a long time.

Clean the appliance with a soft damp cloth.

The appliance will switch off automatically if it is not used for a few seconds and BFM mode is selected.

IMPORTANT

- This appliance is manufactured in compliance with current recycling regulations.
- Do not dispose of protective packaging material with household waste, but take it to a collection facility.
- Use only original Polti spare parts.

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ENGLISH

Polti products are guaranteed to the original purchaser for two years from the date of purchase as shown on a valid receipt issued by the seller.

This guarantee does not affect the consumer's rights under European Directive 99/44/CE regarding certain aspects of the sale and guarantee of consumer goods; consumers should address claims under these provisions to the retailer where they purchased the goods.

This guarantee is effective in countries where European Directive 99/44/CE has been implemented. Local guarantee legislation applies in all other countries.

During the guarantee period, Polti will repair products with defects in workmanship or materials at no charge to the customer for either labour or materials.

To have work performed under this guarantee, consumers must contact an Authorised Polti Technical Assistance Centre. A list of these

centres and their addresses is attached. The consumer must present a valid receipt issued by the seller at the time of purchase as evidence of the date of purchase.

The guarantee does not cover defects due to incorrect use or use other than that identified in the instruction booklet, which is an integral part of the contract of sale for the product. The guarantee does not cover defects resulting from chance (fire, short circuit) or attributable to a third party (tampering with the product).

The guarantee does not cover parts that are damaged by wear or subject to normal wear (consumables such as filters, brushes, etc.).

Polti will accept no responsibility for any injury to people or animals or damage to objects directly or indirectly resulting from failure to follow the instructions provided in the instruction booklet, in particular the warnings and product maintenance instructions.