
Panasonic

Operating Instructions

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Sonic Vibration Toothbrush
Schallzahnbürste

Model No. EW1031

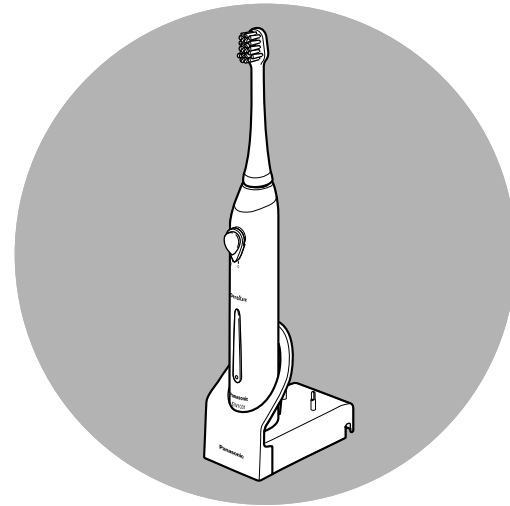
Panasonic Electric Works,(Thailand) Ltd.
Navanakorn Industrial Estate, Zone3 Pathum Thani 12120, Thailand


Before operating this device, please read these instructions completely
and save this manual for future use.

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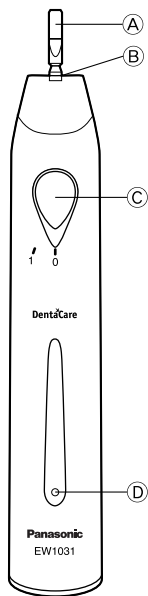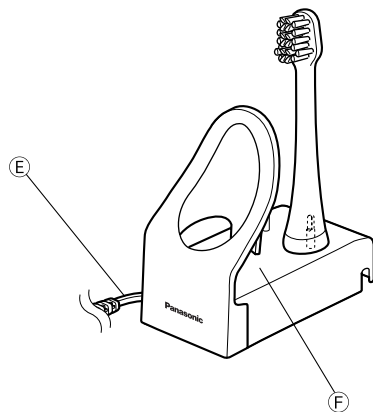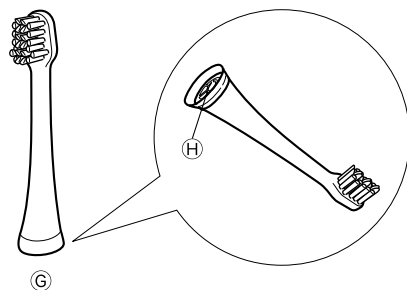
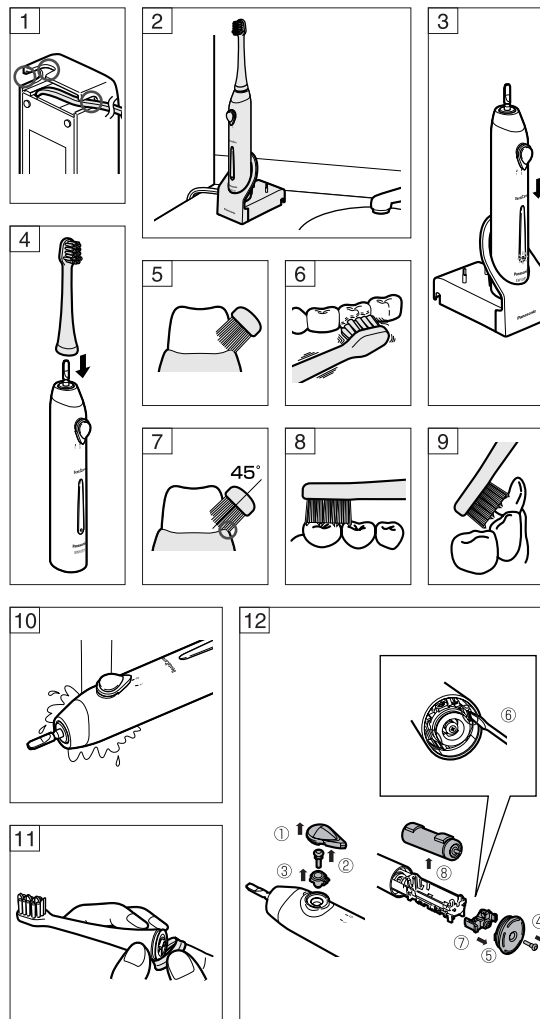
For your personal dental care

Thank you for purchasing the Sonic Vibration Toothbrush.

The Sonic Vibration Toothbrush is a new kind of electric toothbrush designed for those who wish to keep their teeth clean and healthy. This product will enable you to maintain good dental hygiene, which is as important as other forms of personal hygiene.

Sonic Vibration Toothbrush Technology

The Sonic Vibration Toothbrush incorporates innovative motor technology that uses magnetic attraction and repulsion to provide outstanding performance. The brush's 26,000 strokes per minute help remove plaque gently yet effectively.

The intensity of the toothbrush vibration is automatically weakened when the brush is pressed too hard against your teeth.

PARTS IDENTIFICATION

- (A) Toothbrush shaft
- (B) Waterproof rubber seal
- (C) On/off switch
- (D) Charging indicator
- (E) Power cord
- (F) Charger
- (G) V-shaped brush
- (H) Identification ring

BEFORE USE

This high-speed, vibrating, electric toothbrush is different from conventional electric toothbrushes. Therefore:

- You only need to lightly press the bristles against your teeth for brushing.
- When first using this product, some may experience a tickling or tingling sensation caused by the high-speed vibrations. The sensation should subside after you use the toothbrush several more times.

Installing the charger

Cord (Fig. 1)

Place the cord in either the right or left cord guide, depending on where the charger will be placed.

Installing near a sink or on a shelf

Place the charger on a level, secure surface. Insert the plug into an electrical outlet. (Fig. 2)

CAUTION:

- Do not use the charger in a shower stall or other area with running water that come in direct contact with the charger.
- Do not remove the plug from the electrical outlet with wet hands.
- Be sure to use only a 100-240 VAC power supply.

CHARGING THE UNIT

Be sure to charge the main unit for at least 15 hours:

- 1) before using for the first time;
- 2) when you have not used the toothbrush for more than three months; or
- 3) if the toothbrush runs out of power during use.

Be sure to turn off the toothbrush before charging.

Place the toothbrush upright in the charger. (Fig. 3)

The charging indicator will light up.

- The toothbrush charges whenever it is placed in the charger.
- The indicator light does not go off even after the toothbrush is fully charged.
- The toothbrush will provide about 30 minutes of use when the battery is fully charged.
- The operating duration will become shorter as the battery gets older.

CAUTION:

Do not leave metal objects, such as coins or paper clips, on top of the charger. They may become very warm.

NOTE:

- During charging, the toothbrush will get warm, but this is normal.
- Use only the exclusive charger provided.
- If the toothbrush operates only for a short time after being fully charged, this means that the battery is nearing the end of its life cycle. Take the unit to an authorised service centre to have the battery replaced.

HOW TO USE

1. Place the brush head firmly onto the toothbrush shaft. (Fig. 4)
 - There is a small space between the brush and the main unit. This is normal and is necessary for the head to vibrate properly.
2. Place the toothbrush in your mouth.
 - Apply small amount of normal toothpaste to the toothbrush, if desired, before turning the power on to prevent the toothpaste from spattering.
3. Turn the on/off switch to "1" position to turn on the toothbrush.
4. When you have finished brushing, turn the on/off switch to "0" position to turn off the power.
 - To prevent foam from spattering, turn the toothbrush off before removing it from your mouth.

Strong pressure preventive function

Pressing the brush bristles too hard against your teeth will activate this feature, causing the intensity of the vibration to be weakened automatically as a safeguard against excessive pressure.

PRECAUTION:

You may experience slight bleeding of the gums for the first 2 or 3 days of using this toothbrush, but such bleeding should stop. If the bleeding persists for more than a week, stop using this toothbrush and consult your dentist.

NOTE:

Replacing the brush head every 3 months is recommended.

If the bristles are bent to one side, the toothbrush is not very effective for brushing. Replace the brush with a new one when the bristles are bent or worn.

TIPS ON USAGE

V-Shaped Brush

Brush carefully between teeth or around teeth that are not straight. (Fig. 6)

Place the bristles of the toothbrush gently against your teeth.

Front teeth, and between teeth and gums:

- Place the bristles at 45°. (Fig. 7)

Brush the chewing surface of teeth:

- Brush at a 90° angle. (Fig. 8)

Brush the back of teeth or where

teeth are not straight:

- Brush vertically.
- Move the toothbrush to match the angle of the teeth. (Fig. 9)

AFTER USE

Toothbrush

1. Wash the toothbrush with water and shake well to remove the water. (Fig. 10)
2. Return the brush to the brush stand.
3. Use a cloth to clean the main unit. Wash hard-to-remove stains with water.
 - Do not submerge the main unit in water for a long time.
 - After washing with water, wipe dry with a cloth.
 - Do not scratch or damage the waterproof rubber seal. Doing so may allow water to seep inside and damage the main unit.
4. Place the main unit on the charger.

If the identification ring comes off, put it back on. (Fig. 11)

- Do not forcibly pull the identification ring since it can break easily.
- Using the toothbrush without the identification ring has no effect on performance; the ring is simply for user identification purposes.

MAINTENANCE

Charger

Use a cloth to wipe off the charger.

- If cleaning requires a detergent, use a mild, neutral cleanser.

- Do not use harsh chemicals such as thinner, benzene or alcohol as these will discolour the surface.

About once every six months wipe the prongs of the power plug with a dry cloth.

CAUTION:

- Do not wash or submerge the charger in water.
- Clean the power plug periodically. Dirt and dust buildup may cause electrical shock or reduce the effectiveness of the insulation, thus resulting in overheating and a fire.
- The supply cord can not be replaced. If the cord is damaged, the appliance should be scrapped.

TROUBLESHOOTING

Problem	Possible cause	Remedy
If the toothbrush is pressed hard against your teeth, the intensity of the vibration weakens.	Excessive pressure preventive function is activated.	Place the brush gently against your teeth.
The toothbrush operates only for a short period of time even immediately after being charged.	The main unit was not positioned straight up in the charger and was unable to charge properly. (See page 5.)	Stand the main unit up straight so that the bottom touches the charging base, and make sure the charging indicator comes on.
	The battery has reached the end of its life cycle (about 3 years).	Take the unit to an authorised service centre to have the battery replaced. (Do not attempt to replace the battery yourself since you may damage the unit's waterproofing or internal mechanisms.)
	Charging time was insufficient.	Charge the main unit for at least 15 hours. (See page 5.)
The toothbrush doesn't operate.	You have just purchased the product or it has not been used for more than 3 months.	Charge the main unit for at least 15 hours. (See page 5.)
The direction of the brush movement changes. (The brush head moves from side to side.)		Not a malfunction. Remove the brush head once and put it back to the correct position.

Problem	Possible cause	Remedy
The toothbrush is not working properly. (For example, its movement is uneven.)		Place the main unit in the charger for at least 15 hours, then turn it on again to see if it operates properly.

If the device still does not function properly, stop using it and contact an authorised service centre.

WARNING:

Do not disassemble the product at any time other than when disposing of it. Doing so may result in short circuit, fire, or injury. For repair or maintenance (such as battery replacement), take the product to an authorised service centre. The supply cord of the charging adaptor cannot be replaced; if the cord is damaged, the entire charger should be discarded.

Information on Disposal for Users of Waste Electrical & Electronic Equipment (private households)

This symbol on the products and/or accompanying documents means that used electrical and electronic products should not be mixed with general household waste.

For proper treatment, recovery and recycling, please take these products to designated collection points, where they will be accepted on a free of charge basis.

Alternatively, in some countries you may be able to return your products to your local retailer upon the purchase of an equivalent new product.

Disposing of this product correctly will help to save valuable resources and prevent any potential negative effects on human health and the environment which could otherwise arise from inappropriate waste handling. Please contact your local authority for further details of your nearest designated collection point. Penalties may be applicable for incorrect disposal of this waste, in accordance with national legislation.

For business users in the European Union

If you wish to discard electrical and electronic equipment, please contact your dealer or supplier for further information.

Information on Disposal in other Countries outside the European Union

This symbol is only valid in the European Union.

If you wish to discard this product, please contact your local authorities or dealer and ask for the correct method of disposal.

DISPOSAL OF THE PRODUCT (Fig. 12)

The power source of this product is nickel-metal hydride battery which can be recycled. Please contact an authorised service centre for information on how to recycle this battery.

WARNING:

Do not disassemble the product at any time other than when disposing of it. Failure to observe this point may result in short circuit, fire, or injury.

For repair or maintenance (such as battery replacement), take the product to an authorised service centre.

Removing the Ni-MH battery before disposing of the product

- Before removing the battery, be sure to completely drain the unit of power.
- Use a screwdriver to loosen the screw on the bottom of the product.
- While holding the toothbrush shaft, remove the screw then remove the parts indicated by numbers ①, ②, and ③ in the diagram.
- Be careful not to short-circuit the negative and positive poles of the battery.

Please dispose of the device and the battery in accordance with national or public regulations or environmental protection or recycling regulations.

SPECIFICATIONS

Rated voltage:

100 - 240 V AC

50 - 60 Hz

Battery:

1 Ni-MH battery (1.2 V DC; 1,100 mAh)

Charging time: 15 hours

Input current: 35 mA

Power consumption:

2 W (when fully charged)

Brush rotating speed:

approx. 26,000 strokes/min. (with no load)

Brush movement width:

approx. 0.8 mm

Operating time:

approx. 40 minutes of continuous use (with 300 g load at a room temperature of 20°C after the battery has been fully charged.)

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