KH4011

- (D) (A)** Massagematte
Bedienungsanleitung
- (F)** Tapis d'automassage'
Mode d'emploi
- (I)** Tappetino per massaggi
Istruzioni per l'uso
- (E)** Estera de masaje
Instrucciones de servicio
- (GB) (IRL)** Massage mat
Operating Instructions
- (P)** Esteira massajadora
Manual de instruções
- (NL)** Massagemat
Gebruiksaanwijzing
- (GR)** Στρώμα μασάζ
Οδηγίες χρήσης
- (PL)** Poduszka masażowa
Instrukcja obsługi
- (FIN)** Hierontamatto
Käyttöohje
- (CZ)** Masážni poduška
Návod k obsluze
- (S)** Elektrisk massagedyna
Bruksanvisning
- (SK)** Masážna poduška
Návod na obsluhu
- (N)** Massasjematte
Bruksanvisning
- (H)** Massázspárna
Használati utasítás

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⚠ Notes on Safety ⚠

Keep this manual for future reference
– and pass it on with the appliance to any future users!

To avoid life-threatening electric shocks:

- **Only allow customer service to replace damaged cables, before you continue using the appliance.**
- **Use the appliance in dry areas only.**
Do not use outdoors, in the bathtub or in the shower.
- **Do not touch the massage mat, if it is wet.** Pull out the plug from the power socket.

To prevent the risk of fire or injury:

- **Plug the massage mat in only into a 230 V power socket or a 12 V cigarette lighter.** For this purpose, only use the appropriate adapter supplied with the appliance.
- **Do not allow children and elderly persons to handle the appliance unsupervised.** They are not always able to recognize potential dangers correctly.
- **Consult your doctor before using the massage mat,** in case you are being treated for some ailment.
- **Do not use the heating function if you have a sensitive skin or poor circulation.**

Accident hazard!

- **Operate the massage mat in the car only when the vehicle has been stopped.**

1. Intended Use

This massage mat is intended exclusively for private, household use. It is intended for the relaxation massage of shoulders, back, hips and thighs.

This massage mat is not suitable for any other use. It is not intended for use in medical/ therapeutic or commercial fields.

2. Delivery Contents

- 1 massage mat
- 1 mains adaptor
- 1 12 V adaptor
- 1 Operating Manual
- 1 Warranty card

3. Technical Data

Nominal voltage: 12 V DC / 800 mA

4. Connection

Accident hazard!

- **Operate the massage mat in the car only when the vehicle has been stopped.**

If you want to use the massage mat in the car...

- ⇨ firstly insert the plug of the 12 V adaptors into the remote control jack provided for the purpose...
- ⇨ and then insert the large plug into the cigarette lighter of your car.


If you want to use the massage mat at home...

- ⇒ first of all, insert the small plug of the mains adaptor into the remote control jack provided for the purpose...
- ⇒ and then insert the power plug into the power socket.

5. Operate with the remote control

With the remote control function keys, you can easily control the individual functions of the massage mat:

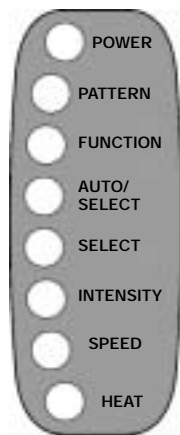PATTERN

You can set the massage pattern with this function key.

When you switch on the massage mat, you receive a **vibrating message**.


Press the PATTERN key...

- ⇒ once for a **tapping message**,
- ⇒ twice for a **kneading message**,
- ⇒ thrice for a **undulating message**.

FUNCTION

You can select the sequence of the massage zones with this function key.

Several zones can be activated at the same time.

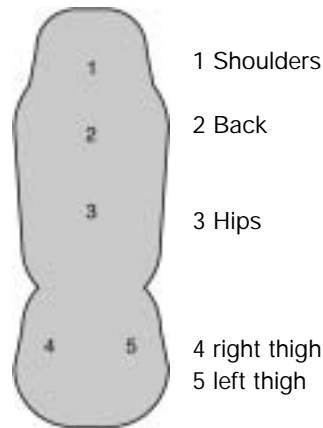POWER

Press this key...

- ⇒ once to switch on the massage mat,
- ⇒ once again to switch off the massage mat.

The massage mat switches off automatically after 15 minutes of operation. To switch on the massage mat again,

- ⇒ press the Power key again.

The zones active at a particular time are displayed with green flashing diodes on the right edge of the remote control. Thus, the topmost diode is for zone 1, the second for zone 2, etc.

When you switch on the massage mat, the standard sequence is as follows:

1 → 2 → 3 → 4 → 5 → 5 → 4 → 3 ...

Press the FUNCTION key...

- ⇒ once for the following sequence:
1+2+3 → 4+5 → 1+2+3 ...

- ⇒ twice for the following sequence:
1+2 → 2+3 → 3+4 → 4+5 →
4+5 → 3+4 ...
- ⇒ thrice for the following sequence:
3 → 2+4 → 1+5 → 1+5 → 2+4 ...

AUTO/SELECT

- ⇒ Keeping this key pressed, interrupts all the massage processes.

SELECT

You can specifically select one or more massage zones with this function key.

- ⇒ During the massage, press the AUTO/SELECT key to retain the zone that is active at the moment.
- ⇒ Press the FUNCTION key to change the massage zones.
- ⇒ Press the FUNCTION key to simultaneously activate several zones.
- ⇒ Press the AUTO/ SELECT key once again to deactivate the function. The massage mat goes back to the standard cycle again.

INTENSITY

You can control the intensity of the massage with this key.

- ⇒ Press this key for the desired massage intensity: high → medium → low

SPEED

You can control the massage motor cycle speed with this key.

- ⇒ Press this key for the desired speed:
Fast → medium → slow

HEAT

- ⇒ Press this key to switch on and switch off the heating function.
- The lower diode on the edge of the remote control flashes red light if the heating function is on.

7. Cleaning

Warning!

- **Pull out the plug from the power socket**, before you clean the massage mat.
- **Never immerse the massage mat in water or other liquids for cleaning.**

- ⇒ Remove fluff with a clothes brush.
- ⇒ You can treat stains with a slightly wet cloth.

9. Disposal

Inquire with your municipal authorities or local government about options for the environment-friendly and appropriate disposal of the appliance.

10. Warranty & Service

The warranty conditions and service address can be found in the enclosed warranty slip.

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