**CTC** **CLATRONIC®**

## **Bedienungsanleitung/Garantie**

Gebbruksaanwyzing/Garantie • Betjeningsvejledning/garanti  
Mode d'emploi/Garantie • Instrucciones de servicio/Garantía  
Manual de instruções/Garantia • Istruzioni per l'uso/Garanzia  
Instruction Manual/Guarantee • Instrukcja obsługi/Gwarancja  
Návod k použití/Záruka • A használati utasítás/garancia  
Руководство по эксплуатации/Гарантия

## **Elektrische Zahnbürste**

Elektrische tandenborstel • Elektrisk tandbørste  
Brosse à dents électrique • Cepillo de dientes eléctrico  
Escova de lavar dentes eléctrica • Spazzolino da denti elettrico  
Electric Toothbrush • Elektryczna szczoteczka do zębów  
Elektrický zubní kartáček • Elektromos fogkefe  
Электрическая зубная щётка

**EZ 2597**

**CE**

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
## General safety instructions

**GB**

- **Before commissioning this device please read the instruction manual carefully.**
  - Please keep these instructions, the guaranty certificate, the sales receipt and, if possible, the carton with the inner packaging!
  - The device is intended exclusively for private and not for commercial use!
  - Always remove the plug from the socket whenever the device is not in use, when attaching accessory parts, cleaning the device or whenever a disturbance occurs. Switch off the device beforehand. Pull on the plug, not on the cable.
  - In order to protect children from the dangers of electrical appliances, never leave them unsupervised with the device. Consequently, when selecting the location for your device, do so in such a way that children do not have access to the device. Take care to ensure that the cable does not hang down.
  - Test the device and the cable regularly for damage. If there is damage of any kind, the device should not be used.
  - Do not repair the device yourself, but rather consult an authorised expert.
  - For safety reasons a broken or damaged mains lead may only be replaced by an equivalent lead from the manufacturer, our customer service department or a similiary qualified person.
  - Keep the device and the cable away from heat, direct sunlight, moisture, sharp edges and suchlike.
  - Never use the device unsupervised! Switch off the device whenever you are not using it, even if this is only for a moment.
  - Use only original accessories.
  - Do not use the device outdoors\*).
  - Under no circumstances must the device be placed in water or other liquid, or come into contact with such\*).
  - Do not use the device with wet or moist hands\*).
  - Should the device become moist or wet, remove the mains plug from the socket immediately\*).
  - Do not reach into the water.
  - Use the device only for the intended purpose.
- \*) Please note the following „Special safety instructions...“.

## Inserting the Batteries

- Turn the cover of the battery compartment located on the bottom of the toothbrush in the direction "Open".
- Pull off the cover and insert a 1.5-volt battery of type AA as shown in the battery compartment.
- Replace the cover on the battery compartment and turn it in the direction "Close".

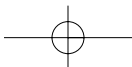**GB**

## Using the Toothbrush

1. Insert a brush on to the handle. The bristles must face in the direction of the switch. It must not be possible to turn the shaft of the brush in this position.
2. Rinse out your mouth thoroughly and wet the bristles. If necessary, you can also use toothpaste. In order to avoid splashing it is better to switch the toothbrush on and off when it is in your mouth.
3. Switch the toothbrush on by pressing the switch to position "I".
4. Move the brush slowly with small circular movements and with slight pressure from tooth to tooth. It is not necessary to press hard. Ensure that the head of the brush also reaches into the spaces between the teeth.
5. After you have finished brushing, rinse out your mouth. Remove the head of the brush and clean it in water.
6. In order to switch off the toothbrush, press the switch to position "0".

Note: It is possible that during the first days of use your gums may bleed. This normally stops after a few days. If the bleeding of the gums does not stop after two weeks, please consult a dentist.

## Cleaning and Care

- Remove the head of the brush from the handle and clean it under running water! Do not use any additives to clean the head of the brush!
- The toothbrush is watertight and can be cleaned under running water! However, it should not be left submerged under water for a prolonged period!
- Do not use any abrasive substances for cleaning! These may damage the seals and the housing!
- Simple soapy water may be used to clean the toothbrush.
- If the bristles lose their stiffness and shape, they should be replaced. Spare toothbrushes for your electric toothbrush can be ordered from the telephone number listed under "Service".

This appliance conforms with the CE directives for radio interference suppression and low-voltage safety and has been built to meet current safety requirements. Subject to technical changes without prior notice!

## Warranty

We undertake to guarantee the appliance sold by our company for 24 months as of the date of purchase (cash receipt).

Within the guarantee period we rectify faults to the machine free of charge which result from material or manufacturing faults, either by repairing the machine or, or our discretion, by replacement. Any such repair or replacement associated with the

guarantee shall not cause the guarantee period to be extended, nor be deemed to be the start of a new guarantee period!

Proof for the warranty is provided by the purchase document. The replacement and repair will not be performed free of charge without this proof of purchase.

Should the warranty be valid, please consign the complete appliance in its original packing together with the cash receipt to your local dealer.

Defects to the accessories do not entitle to the change of the appliance. These will be replaced free of charge. In this case do not send back the appliance. Just order the faulty accessory! Broken glass parts are replaced at the expense of the owner!

Faults to consumption accessories or wearing parts (e.g. engine graphite, saw blades, kneading hooks, drive belts etc), as well as cleaning, maintenance or the replacement of wearing parts, do not fall under the guarantee and are therefore subject to a charge!

The warranty is terminated in case of handling on the part of Third parties.

### **After the warranty**

Once the warranty has expired, repairs can be carried out by specialized technicians or repair services against the payment of the service.

### **Service**

Dear customer,

Wearing parts which are not replaced by a specialist, replacement cans, vacuum cleaner filter bags etc., can be ordered on the following telephone hotline by stating the machine type:

**+49/2152/2006-888.**

For orders by Internet we have set up the following e-mail address:

**clatronic\_hotline@t-online.de**

Are you interested in other Clatronic products?

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